

Combined Follow-up Report

Date of Birth N/R N/R N/R

Device ACCOLADE MRI L311/269927

Last Office Interrogation

27 Mar 2017 Implant Date N/R

My Alerts

There are no alerts to display.

Events Since Last Reset (23 Mar 2017)

23 Apr 2017 19:17 RA Auto 23 Apr 2017 19:16 RV Auto

Battery OK

Approximate time to explant: 9.5 years
Magnet Rate 100 min⁻¹

One Year Remaining

Explant

		Explain =		
Leads Data	Implant N/R N/R N/R	Previous Session	Most Recent	
Atrial				
Intrinsic Amplitude	N/R mV	2.2 mV	2.0 mV	
Pace Impedance	N/R Ω	615 Ω	619 Ω	
Pace Threshold	N/R V @ N/R ms	Auto 1.7 V @ 0.4 ms	Auto 0.9 V @ 0.4 ms	
Ventricular				
Intrinsic Amplitude	N/R mV	14.8 mV	16.9 mV @ 50 min ⁻¹	
Pace Impedance	N/R Ω	687 Ω	733 Ω	
Pace Threshold	N/R V @ N/R ms	Auto 0.7 V @ 0.4 ms	Auto 0.5 V @ 0.4 ms	
Cattings				

Settings	5
----------	---

Ventricular Tachy Settings		Atrial Tachy Settings
Vantainulas Tanku FOM Ctarana	0	ATD Mada Cuitala 4

Ventricular Tachy EGM Storage On ATR Mode Switch 170 min⁻¹ DDI Detection Rate 160 min⁻¹

DDDR

Brady Settings

Mode

RYTHMIQ™	AAIR With	Atrial	Auto 2.0 V @ 0.4 ms
	VVI Backup	Ventricular	Trend 2.5 V @ 0.4 ms
Lower Rate Limit	55 min ⁻¹	Sensitivity	
Maximum Tracking Rate	130 min ⁻¹	Atrial	Fixed 0.5 mV
Maximum Sensor Rate	130 min ⁻¹	Ventricular	Fixed 2.0 mV
Paced AV Delay	80 - 180 ms	Leads Configuration (Pace/So	ense)
Sensed AV Delay	65 - 150 ms	Atrial	Bipolar
A-Refractory (PVARP)	240 - 280 ms	Ventricular	Bipolar
V-Refractory (VRP)	230 - 250 ms		

Rate Adaptive Pacing

Pacing Output

Minute Ventilation On Accelerometer Passive

2869 Software Version: 2.02 L311 Firmware Version: F_v1.00.00 © 2014

Boston Scientific Corporation
or its affiliates. All rights reserved.

Page 1 of 7

Clinician Signature:


Trends					
Today's Activity:	Program and	Interrogate			
R = Counter Reset		Jan 17 ₁	Feb 17	Mar 17	Apr 17
P = Program				P	
I = Interrogate				I I	
F		Jan 17	Feb 17	Mar 17	Ann 17
Events	А		LED 17	Wat 17	Apr 17
	А				
	V				
	Other				
Activity Level		Jan 17	Feb 17	Mar 17	Apr 17
•	20 15				, , ,
	% 10				\
	5	1			V V V
	0				
Heart Rate		Jan 17	Feb 17	Mar 17	Apr 17
Max					
Mean	min ⁻¹				
– Min					
05.1111			F-L 49		
SDANN		Jan 17	Feb 17	Mar 17	Apr 17
	ms				
HRV Footprint		Jan 17	Feb 17	Mar 17	Apr 17
The Footprine					·
	%				
	70				
Autonomic		Jan 17	Feb 17	Mar 17	Apr 17
Balance Monitor					
AT/AF Burden	4	Jan 17	Feb 17	Mar 17	Apr 17
	3	1 ;			
	h 2				
	1				
		Jan 17	Feb 17	Mar 17	Apr 17
V Rate during		14 IIA	LEN T\	widi 1/	Apr 17
AT/AF		VT			
Max	min ⁻¹				
Mean					


2869 Software Version: 2.02 L311 Firmware Version: F_v1.00.00 © 2014

Boston Scientific Corporation
or its affiliates. All rights reserved.

Page 2 of 7

Clinician Signature:	

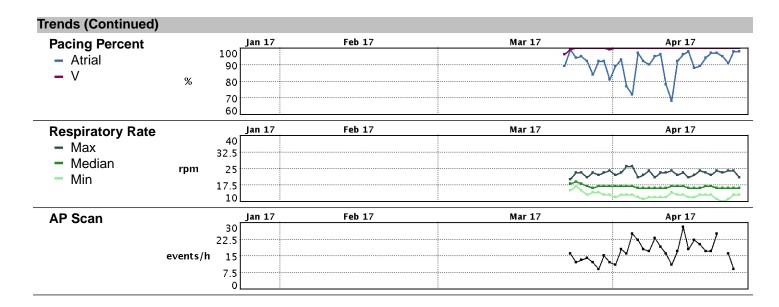

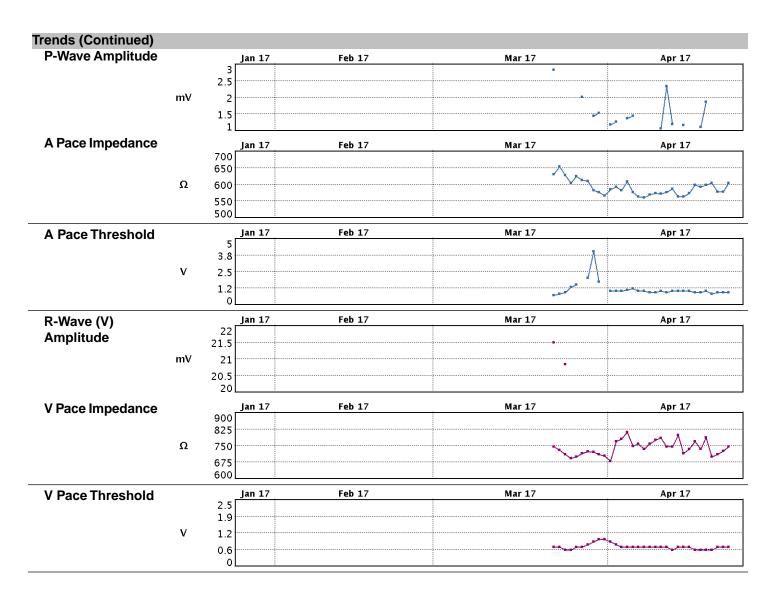

Ventricular Tachy Counters	Since Last Reset 23 Mar 2017	Device Totals
Ventricular Episode Counters		
Total Episodes	0	0
VT Episodes (V>A)	0	0
SVT Episodes (V<=A)	0	0
Nonsustained Episodes	0	0
Brady Counters	Reset Before Last N/R day(s) N/R to 23 Mar 2017	Since Last Reset 32 day(s) 23 Mar 2017 to Today
Counters		
% A Paced	0	91
% V Paced	0	100
Intrinsic Promotion		
AV Search +		
% Successful	0	0
Rate Hysteresis		
% Successful	0	0
Atrial Arrhythmia		
% AT/AF	N/R	0
Total Time in AT/AF (min)	N/R	0.0
Episodes by Duration		
< 1 minute	0	0
1 min - < 1 hr	0	0
1 hr - < 24 hr	0	0
24 hr - < 48 hr	0	0
> 48 hr	0	0
Total PACs	0	1.1K
Ventricular Counters		
Total PVCs	0	1.3K
Three or More PVCs	0	3
AT/AF Overview: Since Last Reset (23 Mar 20	17 to Today)	

AT/AF Overview: Since Last Reset (23 Mar 2017 to Today)

AT/AF: 0% AT/AF Events: 0 Total Time in AT/AF (min): 0.0 ATR Episodes initiated within 2 hours of lead attach are not stored in Arrhythmia Logbook.


AT/AF Overview: Reset Before Last (N/R to 23 Mar 2017)


ATR Episodes initiated within 2 hours of lead attach are not stored in Arrhythmia Logbook.


2869 Software Version: 2.02 L311 Firmware Version: F_v1.00.00

Clinician Signature:


Histograms Reset Before Last (N/R)

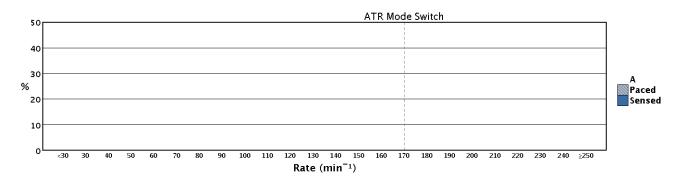

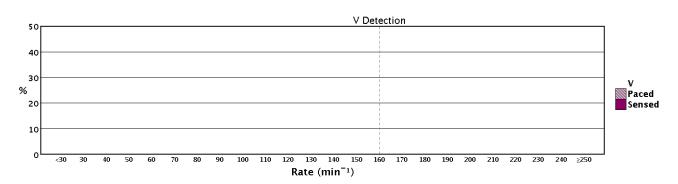

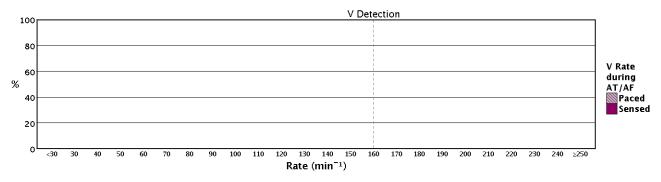

Histograms Since Last Reset (23 Mar 2017)

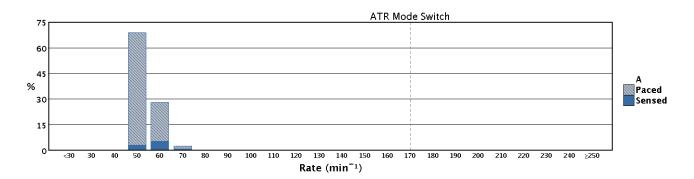


2869 Software Version: 2.02 L311 Firmware Version: F_v1.00.00 © 2014


Boston Scientific Corporation
or its affiliates. All rights reserved.

Page 6 of 7

Clinician Signature:		

Histograms Since Last Reset (23 Mar 2017) (Continued)

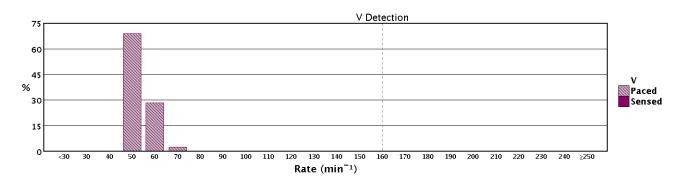

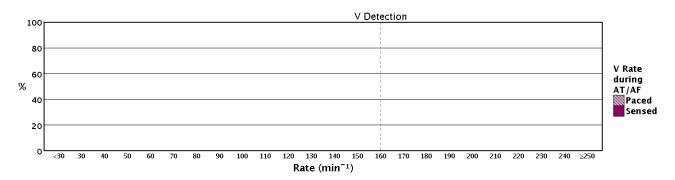

Heart Rate Variability

Last Measured:

HRV data unavailable because device was reprogrammed in the last 24 hours.

Reference:

HRV data unavailable because device was reprogrammed in the last 24 hours.