

FastPath™ Summary

 **2 Alerts**

Page 1 of 1

Battery


Longevity: 5,1-5,4 yrs



Implant Date:	8 Mar 2022
Last Max Charge	8,6 sec (Therapy)
Battery Current	21 uA
Remaining Capacity to ERI	91%

Test Results 9 Mar 2022



 Automatic

	Capture	Sense	Lead Impedance
A	Not Performed 0,75V @ 0,7ms (Bi)	Not Performed None (Bi)	Not Performed 390 Ω (Bi)
RV	Not Performed 0,5V @ 0,5ms (Bi)	Not Performed 6,0mV (RV Bi)	Not Performed 410 Ω (Bi)
LV	Not Performed 1,0V @ 0,5ms (D1-M2)		Not Performed 750 Ω (D1-M2)
HV			40 Ω (RV to SVC & Can)  37 Ω (RV to SVC & Can)

Parameters

Mode	DDDR	Zone Configuration	VT	VF
Base Rate	▶ 90 bpm	Detection Criteria	171 bpm	214 bpm
Max Track Rate	110 bpm	Therapy (ENABLED)	ATP x3	ATP x1
Paced AV Delay	200 ms		25,0 J	36,0 J
Sensed AV Delay	150 ms		36,0 J	40,0 J
Ventricular Pacing	LV → RV, 30 ms		40,0 J x2	40,0 J x4

Capture & Sense

	A	RV	LV	LV2
Cap Confirm	Monitor	Off	Off	
Pulse Amplitude (Margin)	2,5 V (1.3:1)	2,5 V	2,5 V	
Pulse Width	0,5 ms	0,5 ms	0,5 ms	
AutoSense	On	On		
Sensitivity	Auto 	Auto 		

Diagnostics Summary

Since 9 Mar 2022

AP	98 %
BP	98 %
AMS Episodes	7
Mode Switch	<1%
AT/AF Burden	<1%

VT/VF Episodes: 3

Since 9 Mar 2022

	VT	VF
Episodes	0	3
ATP Delivered	0	3
Shocks Delivered	0	4
SVT Episodes: 0		
Non-sustained Episodes: 0		

Alerts



Episodes with alert conditions (3)
A. Pulse Amplitude Safety Margin < 2:1

▶ Manual-programmed ↳ Auto-programmed

Episodes Summary

Page 1 of 1

 Episodes Last Cleared 8 Mar 2022 13:36
 SEGMs Last Cleared 8 Mar 2022 13:18

Last Read 9 Mar 2022 13:37

Therapy Summary

	VT	VF
ATP Delivered	0	5
Shocks Delivered	0	6
Max Energy Shocks	0	2


Last HV Lead Impedance 40 Ω

Total Aborted Shocks 1

Results of ATP Delivery

	VT	VF
Episodes Terminated	0	1
Episodes Not Terminated	0	4
Accelerations	0	

Episode Tree

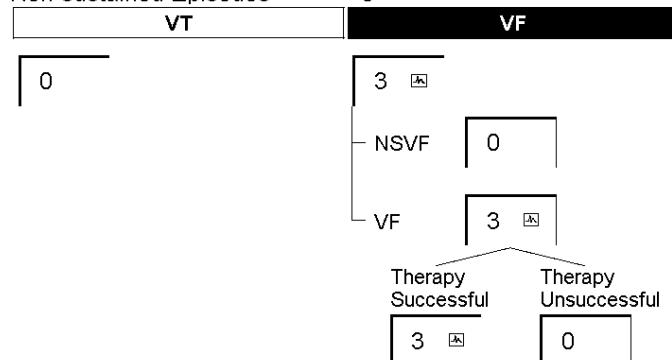 Total VT/VF Episodes 3
 SVT Episodes 0
 Non-sustained Episodes 0


VT/VF Episodes

Date / Time	Type	Rate (bpm)	Duration (M:S)	Therapy Delivered	Alerts
9 Mar 2022 13:21	VF	240	00:29	ATP, 36J, 40J	⚠ x2
9 Mar 2022 13:21	VF	235	00:17	ATP, 36J	⚠ x1
9 Mar 2022 13:20	VF	235	00:17	ATP, 36J	⚠ x1

Other Episodes

No episodes recorded since last session

Extended Diagnostics


Page 1 of 4

Parameters

Mode	DDDR	Max Track Rate	110 bpm
Base Rate	70 bpm	Paced/Sensed AV Delay	200/150 ms
Max Sensor Rate	110 bpm		

Events

Since 9 Mar 2022

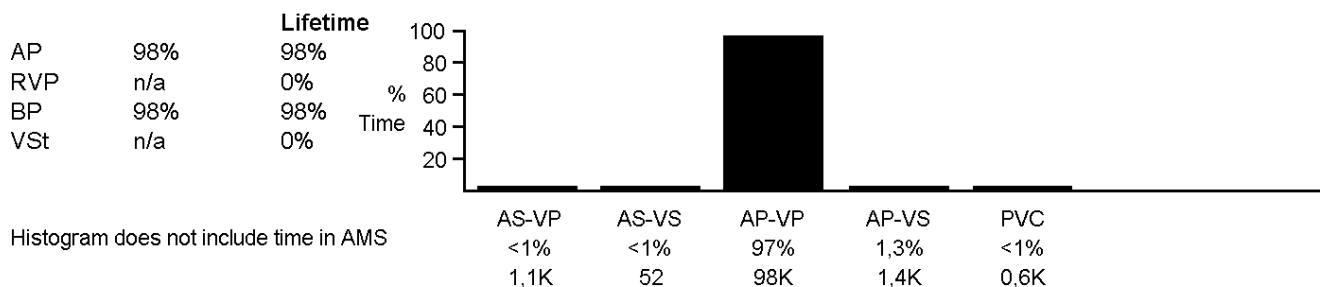

AV Interval

Since 9 Mar 2022

AV Conduction

<1 %


VS-AS Intervals

0

< 80 ms

■ (AS-VP, AP-VP)

□ (AS-VS, AP-VS)

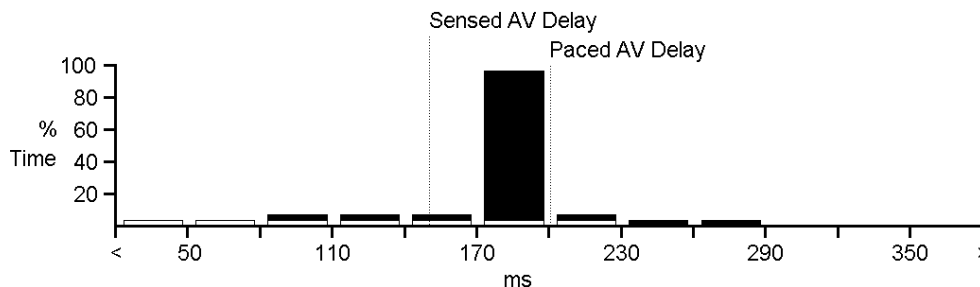

Excludes time in AMS

AS-VP %			<1	<1	<1	<1	<1	<1	<1				
AP-VP %			<1	<1	<1	98	<1	<1	<1				
AS-VS %	<1	<1	<1	<1	<1	<1	<1	<1	<1				
AP-VS %		<1	<1	<1	<1	<1	<1	<1	<1				
Total %	<1	<1	<1	<1	<1	99	<1	<1	<1				

Atrial Heart Rate Histogram

Since 9 Mar 2022

AF Suppression™

n/a

Time at Max Track


<1%

Rate

■ Paced (AP)

□ Sensed (AS)

○ Sensor-Indicated Rate

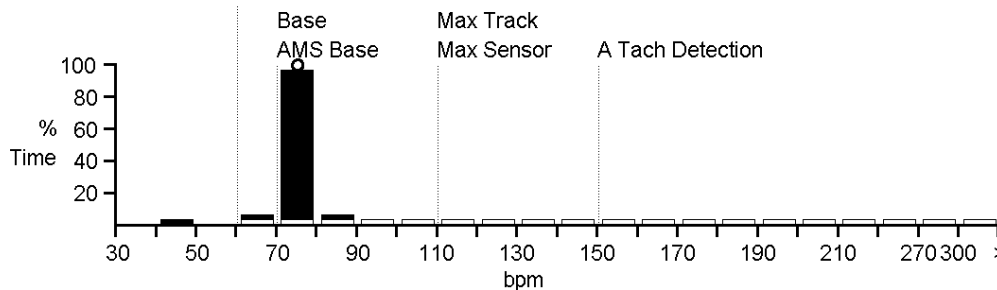

AP %		<1	<1	99	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
AS %		<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Total %		<1	<1	>99	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Sensor %				100															

Page 2 of 4

Since 9 Mar 2022

Rest

Base

AMS Base

Max Sensor

Max Track

%

Time

100

80

60

40

20

30 50 70 90 110 130 150 170 190 210 230 250 270 300 >

bpm

[illegible]

Since 9 Mar 2022

Since 9 Mar 2022

Heart Rate (bpm)	White (%)	Black (%)	Total (%)
70-90	25	55	80
90-110	18	0	18
110-130	5	0	5
130-150	0	0	0
150-170	5	0	5
170-190	5	0	5
190-210	5	0	5

VP %			55							
VS %			26	16	<1			<1	<1	1,1
Total %			81	16	<1			<1	<1	1,1



Extended Diagnostics

Page 3 of 4

AMS Summary

Since 9 Mar 2022

AMS Episodes



7

Mode Switch

<1%

Duration

Peak A Rate

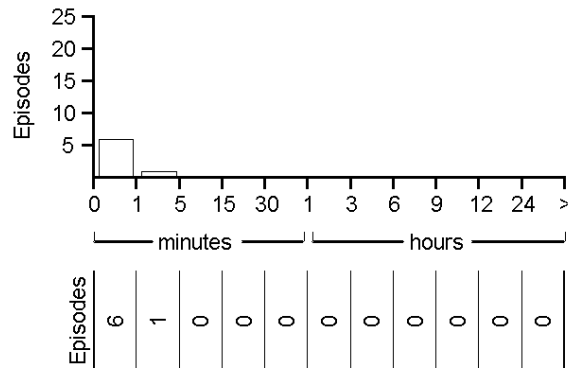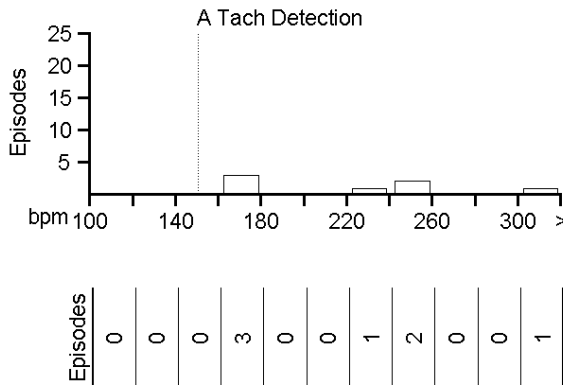

Key AMS Log Episodes	Date and Time	Peak A Rate	Peak V Rate	Duration (D:H:M:S)
Most Recent Episode	9 Mar 2022 7:51	240 bpm	98 bpm	0:00:00:24
Peak V Rate Episode	9 Mar 2022 7:51	240 bpm	98 bpm	0:00:00:24
Longest Episode	8 Mar 2022 13:30	307 bpm	86 bpm	0:00:01:50


AMS Log
EGMs

Date and Time	Peak A Rate	Peak V Rate	Duration (D:H:M:S)
9 Mar 2022 7:51	240 bpm	98 bpm	0:00:00:24
8 Mar 2022 16:32	160 bpm	70 bpm	0:00:00:08
8 Mar 2022 15:57	169 bpm	70 bpm	0:00:00:08
8 Mar 2022 13:33	167 bpm	71 bpm	0:00:00:22
8 Mar 2022 13:33	229 bpm	76 bpm	0:00:00:10
8 Mar 2022 13:32	248 bpm	80 bpm	0:00:00:50
8 Mar 2022 13:30	307 bpm	86 bpm	0:00:01:50

Daily Exercise Training (Activity >1 minute)


Max Heart Rate 143 bpm

- >80% Max HR
- ▨ 60-80% Max HR
- (Target HR)
- 40-60% Max HR

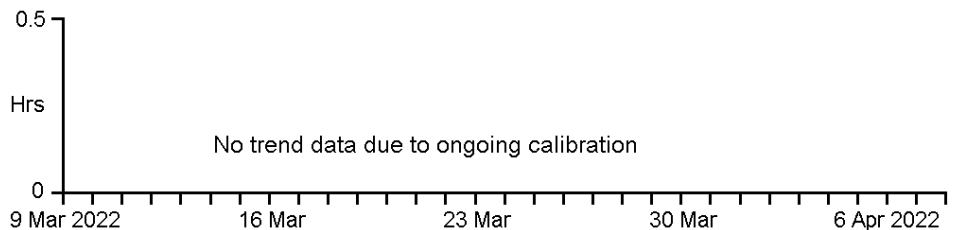

Total Daily Activity

- Activity

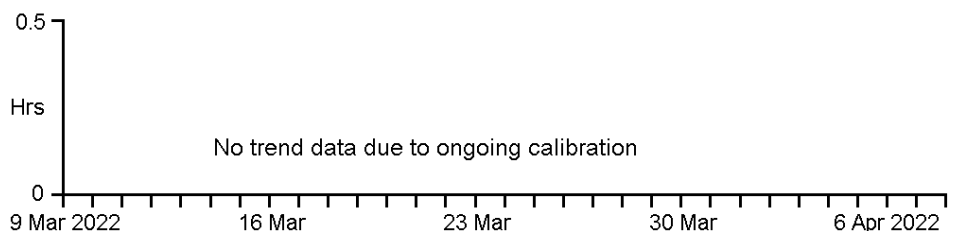

0m 1d 0h 5m Sampled since 8 Mar 2022


0m 0d 0h 3m Sampled in AMS

Last Read: 9 Mar 2022

Extended Diagnostics

Page 4 of 4

CorVue™ Congestion Monitoring is OFF (3 Month View)

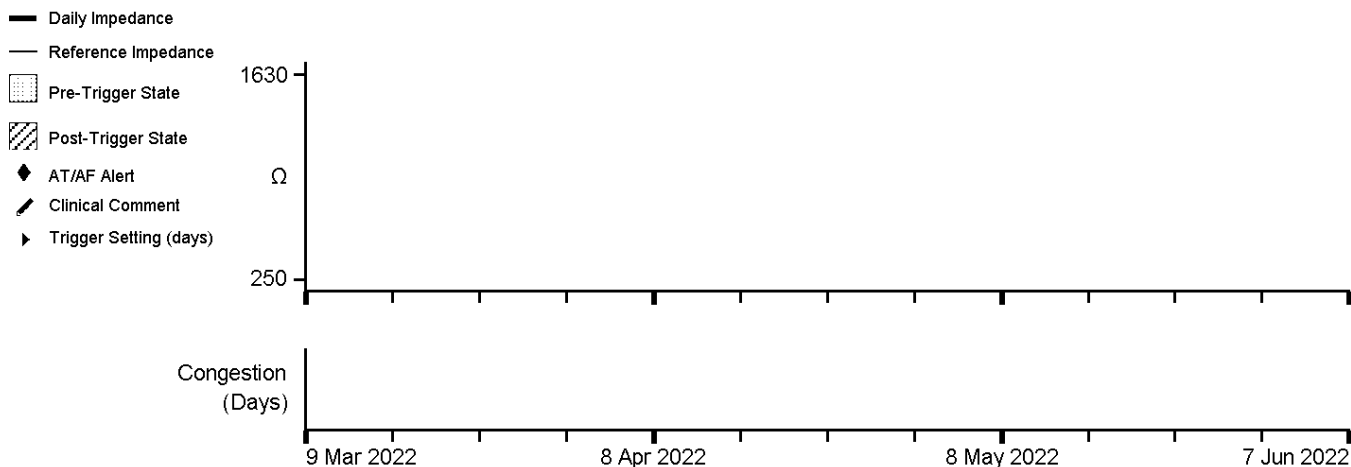

CorVue™ Congestion Monitoring is an additional information source for patient management. It is not meant to replace clinical assessments that are part of standard clinical practice.

Congestion Details (Last Session: 9 Mar 2022(8:37))

No Congestion Episodes recorded

Parameters

Page 1 of 5

Patient

Date of Birth 8 Jun 1944
EF % Unknown

Indications for Implant

Device	Manufacturer	Model	Serial	Implant Date
CRT-D	St. Jude Medical	Quadra Assura MP™ 3371-40QC	9973677	8 Mar 2022
A Lead				
RV Lead	St. Jude Medical	Durata™ 7120Q / 58 cm	BNP078743	8 Mar 2022
LV Lead	St. Jude Medical	Quartet™ 1457Q / 86 cm	DBU030198	8 Mar 2022

Basic Operation

Mode DDDR
Ventricular Pacing LV → RV, 30 ms
V. Triggering Off
Magnet Response Normal
V. Noise Reversion Mode Pacing Off
Episodal Pacing Mode DDI
Sensor On
Threshold (Measured Avg.) Auto (+0.0) (2,0)
Slope 8
Max Sensor Rate 110 bpm
Reaction Time Fast
Recovery Time Medium

Refractories & Blanking

PVARP 275 ms
Post-Vent. Atrial Blanking 90 ms
Rate Responsive PVARP/V Ref Low
Shortest PVARP/V Ref 225 ms
A/V Pace Refractory 190/250 ms
A/V Sense Refractory 93/125 ms
Ventricular Blanking 52 ms
Ventricular Safety Standby On
Arrhythmia Unhiding 3 intervals
PVC Response Off
PMT Response Atrial Pace
PMT Detection Rate 110 bpm

Rates

Base Rate ▶ **90 bpm**
Rest Rate ▶ **Off**
Max Sensor Rate 110 bpm
Max Track Rate 110 bpm
Hysteresis Rate Off
2:1 Block Rate 164 bpm

AT/AF Detection & Response

Auto Mode Switch DDIR
AMS V. Triggering Off
A. Tachycardia Detection Rate 150 bpm
AMS Base Rate ▶ **90 bpm**
AF Suppression™ Off

Delays

Paced AV Delay 200 ms
Sensed AV Delay 150 ms
Rate Responsive AV Delay Medium
Shortest AV Delay 100 ms
SyncAV™ CRT Delta Off

Capture & Sense

	A	RV	LV	LV2
Cap Confirm	Monitor	Off	Off	
Backup Pulse Configuration	Bipolar			
Search Interval	8 hours			
Pulse Amplitude (Margin)	2,5 V (1.3:1)	2,5 V	2,5 V	
Pulse Width	0,5 ms	0,5 ms	0,5 ms	
AutoSense	On	On		
Sensitivity	Auto A	Auto A		

Last Programmed: Today 13:39
Parameters that are "n/a" are not shown

Bold values were changed this session (See
Wrap-up™ Overview report for details)

▶ Manual-programmed **A** Automatic
↳ Auto-programmed

Parameters

Page 2 of 5

Leads	A	RV	LV	LV2
Lead Type	Bipolar	Bipolar	Bipolar	Quadrupolar
Pulse Configuration	Bipolar	Bipolar	Distal tip 1 - Mid 2	
Sense Configuration	Bipolar		RV Bipolar	
Lead Monitoring	Monitor	Monitor	Monitor	
Lower Limit	200 Ω	200 Ω	200 Ω	
Upper Limit	2.000 Ω	2.000 Ω	2.000 Ω	
HVLI Lower Limit		20 Ω		
HVLI Upper Limit		125 Ω		

ShockGuard™ Settings (Zone Configuration)

	VT	VF
Detection Criteria	171 bpm/350 ms 50 intervals	214 bpm/280 ms 12 intervals
SVT Discrimination	On	
Therapy	ATP x3 25,0 J/704 V 36,0 J/845 V 40,0 J/890 V x2	ATP x1 36,0 J/845 V 40,0 J/890 V 40,0 J/890 V x4
VT Therapy Timeout	Off	

SVT Discrimination

SVT Discrimination	Dual Chamber
SVT Discrimination Timeout	Off
SVT Upper Limit	Same as VF

Rate Branch	Additional Discriminators				Diagnosis
AF/A Flutter V < A	Morphology Morphology	Off	Interval Stability Interval Stability Stability Delta AV Association Delta Window Size	On w/AVA 40 ms 60 ms 12 intervals	Deliver therapy only if the criteria in this rate branch indicate VT.
Sinus Tach V = A: AV Interval Delta:	Morphology Morphology	On Off	Arrhythmia Onset Chamber Onset Sudden Onset	On Off	Deliver therapy only if the active criteria in this rate branch indicate VT.

VT/VF V > A	If Ventricular Rate is greater than the Atrial Rate, deliver therapy.
-----------------------	-----------------------------------------------------------------------

Parameters

Page 3 of 5

SecureSense™ Settings

SecureSense™	On
Far Field MD / SecureSense™ Configuration	RV Coil-Can
Timeout Until Therapy	Off
Trigger alert for Non-sustained RV Oversensing (NSO)	5 episodes

ATP Details

ATP Pulse Amplitude	7,5 V
ATP Pulse Width	1,0 ms

	VT		VF
	Therapy 1	Therapy 2	Therapy 1
ATP Type			ATP While Charging
ATP Upper Rate Cutoff			250 bpm/240 ms
Number of Bursts	3		1
Number of Stimuli	8		8
Add Stimuli per Burst	Off		
Burst Cycle Length	85 %		85 %
Min. Burst Cycle Length	200 ms		170 ms
Readaptive	Off		
Scanning	On (Dec)		
Scan Step	10 ms		
Ramp	Off		Off

SenseAbility™ Settings

	A		RV	
			Pacemaker	Defib
Max Sensitivity	0,3 mV		Same as Defib	0,5 mV
	(0,3 mV @ 110 bpm)		(0,5 mV @ 110 bpm)	(0,5 mV @ 110 bpm)
Low Frequency Attenuation				On
	Post-Sensed	Post-Paced	Post-Sensed	Post-Paced
Decay Delay	0 ms	0 ms	60 ms	95 ms
Threshold Start	50 %	0,8 mV	50 %	1,0 mV

DeFT Response™ Settings (Shock Waveform)

Waveform	Biphasic	
Waveform Mode	Tilt	
Shock Configuration	RV to SVC & Can	
RV Polarity	Anode(+)	
DynamicTx™ Over-Current	On	
Detection Algorithm		
	1st Phase	2nd Phase
VF Shocks (Defib)	65 %	65 %
VT Shocks (CVRT)	Same as Defib	Same as Defib
Estimated Defib Pulse Width	5,0 ms	
Estimated CVRT Pulse Width	5,0 ms	
Shock Impedance	40 Ω (Measured)	

 Last Programmed: Today 13:39
 Parameters that are "n/a" are not shown

 Bold values were changed this session (See
 Wrap-up™ Overview report for details)

 ▶ Manual-programmed **A** Automatic
 ↳ Auto-programmed

Parameters

Page 4 of 5

Post-Shock Pacing

Mode	DDD
Base Rate	60 bpm
Pause	2 sec
Duration	0,5 min

Redetection & Post-Detection

VT Redetection	6 intervals
Sinus Redetection	Nominal (5 intervals)
Post VF Detection	Same as VT
Post VT Detection	Same as VT

	A	BV
Pulse Amplitude	7,5 V	7,5 V
Pulse Width	1,5 ms	1,5 ms

Capacitor Maintenance

Charge Interval	4 months
-----------------	----------

MRI

MRI Mode	DOO
MRI Base Rate	85 bpm
MRI Paced AV Delay	120 ms
MRI Atrial Pulse Amplitude	5,0 V
MRI RV Pulse Amplitude	5,0 V
MRI Atrial Pulse Width	1,0 ms
MRI RV Pulse Width	1,0 ms
MRI Atrial Pulse Configuration	Bipolar
MRI RV Pulse Configuration	Bipolar
HV Lead Type	Dual Coil
MRI V. Pacing Chamber	RV Only

Congestion Monitoring

Congestion Monitoring	Off
-----------------------	-----

Trigger Alerts When

	Show on FastPath™	Notify Patient
AT/AF Episode (3 hours)	On	Off
AT/AF Burden (6 hours evaluated daily)	On	Off
V Rate during AT/AF (100 bpm for 6 hours, evaluated daily)	On	Off
A Lead Impedance Out of Range (200 - 2.000 Ω)	On	On
RV Lead Impedance Out of Range (200 - 2.000 Ω)	On	On
Therapy inhibits due to RV lead noise	On	On
5 episodes of Non-sustained RV Oversensing detected	On	Off
LV Lead Impedance Out of Range (200 - 2.000 Ω)	On	On
HV Lead Impedance Out of Range (20 - 125 Ω)	On	On
Percent BiV Pacing Alert (<90% over 7 days)	On	Off
Device at ERI	On	On
Charge Time Limit Reached (32 sec)	On	On
Possible HV Circuit Damage	On	On
Device Parameter Reset	On	On
Backup VVI	On	On

Patient Notifier Delivery

Vibration Duration	6 sec
Number of Notifications	4
Time Between Notifications	10 hours

Parameters

Page 5 of 5

Episode Triggers

Atrial Episode (AMS)	Low
VT (Detection)	High
VF (Detection)	High
Non-sustained RV Oversensing	High
Non-sustained VT/VF	Low
PMT	Off
Noise Reversion	Off
Magnet Response	Low

Stored EGM Configuration

VT/VF EGM Max Duration	1 min
VT/VF Pre-Trigger Max Duration	14 sec
Channels (Storage)	2 (22 min)
Channel 1	A Sense Amp
Channel 2	V Sense Amp

**Episode: VF (240 bpm / 250 ms)**

VT/VF Episode 6 of 6

Page 1 of 5

9 Mar 2022 13:21

Duration 00:29 (M:S)

Detection Criteria 214 bpm

Alerts



At least one shock unsuccessful

Three or more VT/VF episodes in 24 hours

Therapy

- 1) ATP x 1
- 2) Defib 36,0 J (844V)
- 3) Defib 40,0 J (890V)

Results

- 1) VF
- 2) VF
- 3) Below Rate Detection (CL > 995 ms)

HV Therapy

Last HV Lead Impedance	40 Ω
First Charge Time	6,9 sec
Last Charge Time	8,6 sec
Delivered PW	+5,0 ms, -5,0 ms

ATP Therapy Details**VF ATP Therapy****Therapy**

Burst 1 216, 216, 216, 216, 216, 216 ms

Diagnosis Summary

Diagnosis

VF

Time to Diagnosis 2,75 sec

Rate (CL) 240 bpm (250 ms)

DeFT Response™ Settings (Shock Waveform)

			1st Phase	2nd Phase
Waveform	Biphasic	VF Shocks (Defib)	65 %	65 %
Waveform Mode	Tilt	VT Shocks (CVRT)	Same as Defib	Same as Defib
Shock Configuration	RV to SVC & Can			
RV Polarity	Anode(+)			
Defib Max Sensitivity	0.5 mV			

Episode: VF (240 bpm / 250 ms)

VT/VF Episode 6 of 6

Page 2 of 5


9 Mar 2022 13:21



1: A Sense Amp AutoGain (3,4 mm/mV)
 2: V Sense Amp AutoGain (0,7 mm/mV)
 3: Discrimination AutoGain (2,6 mm/mV)

4: Markers

Sweep Speed: 25 mm/s

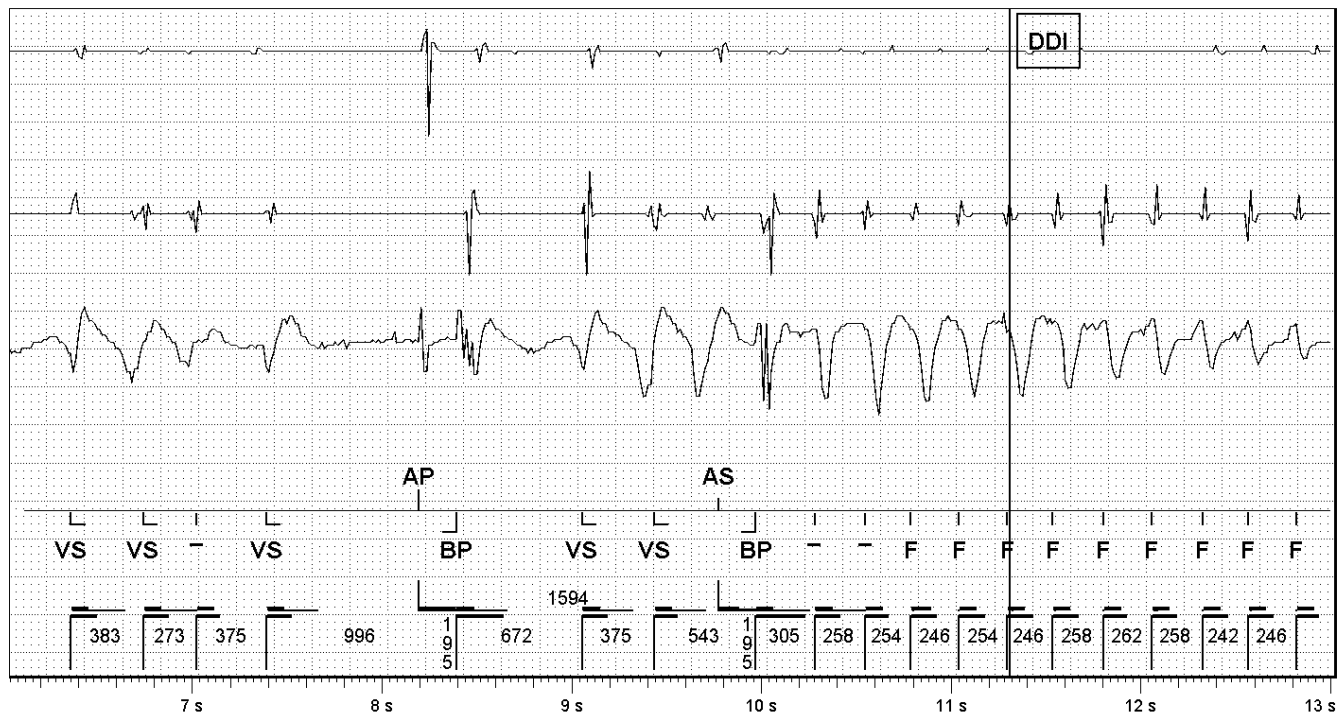




Episode: VF (240 bpm / 250 ms) (Continued)

VT/VF Episode 6 of 6

Page 3 of 5

9 Mar 2022 13:21

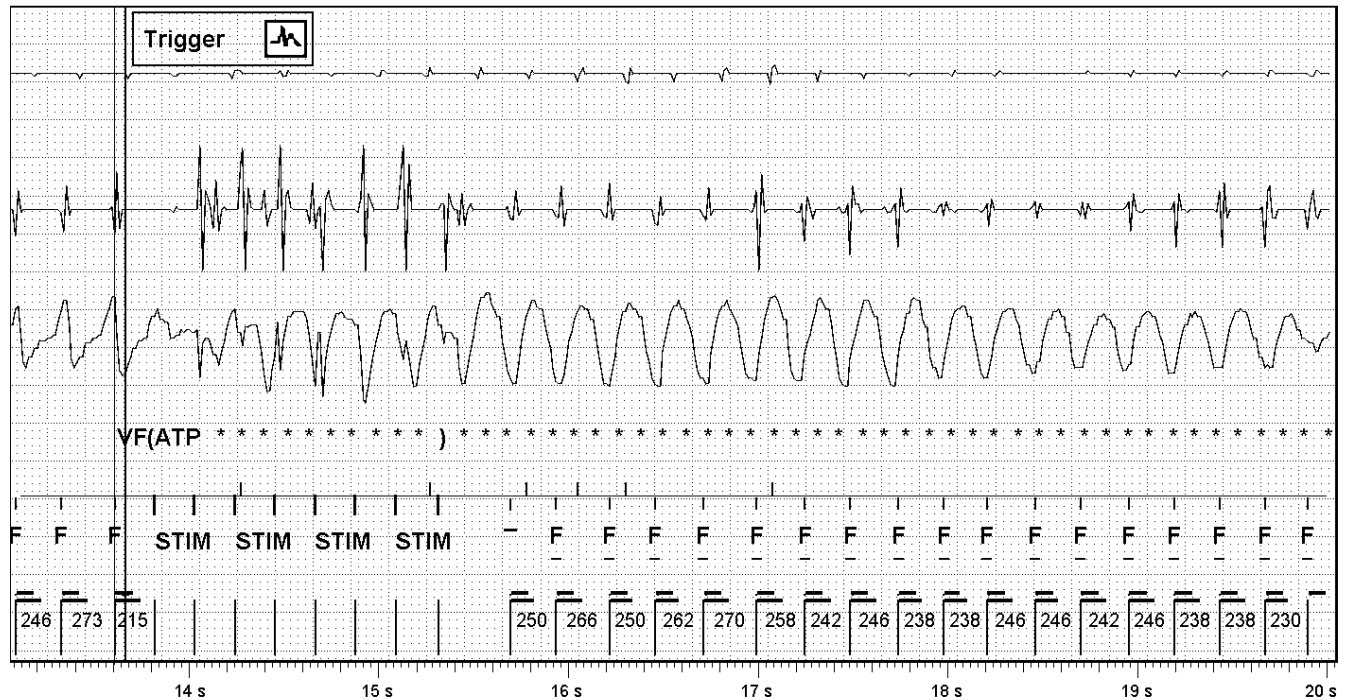


1: A Sense Amp AutoGain (3,4 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,7 mm/mV)

3: Discrimination AutoGain (2,6 mm/mV)

Sweep Speed: 25 mm/s

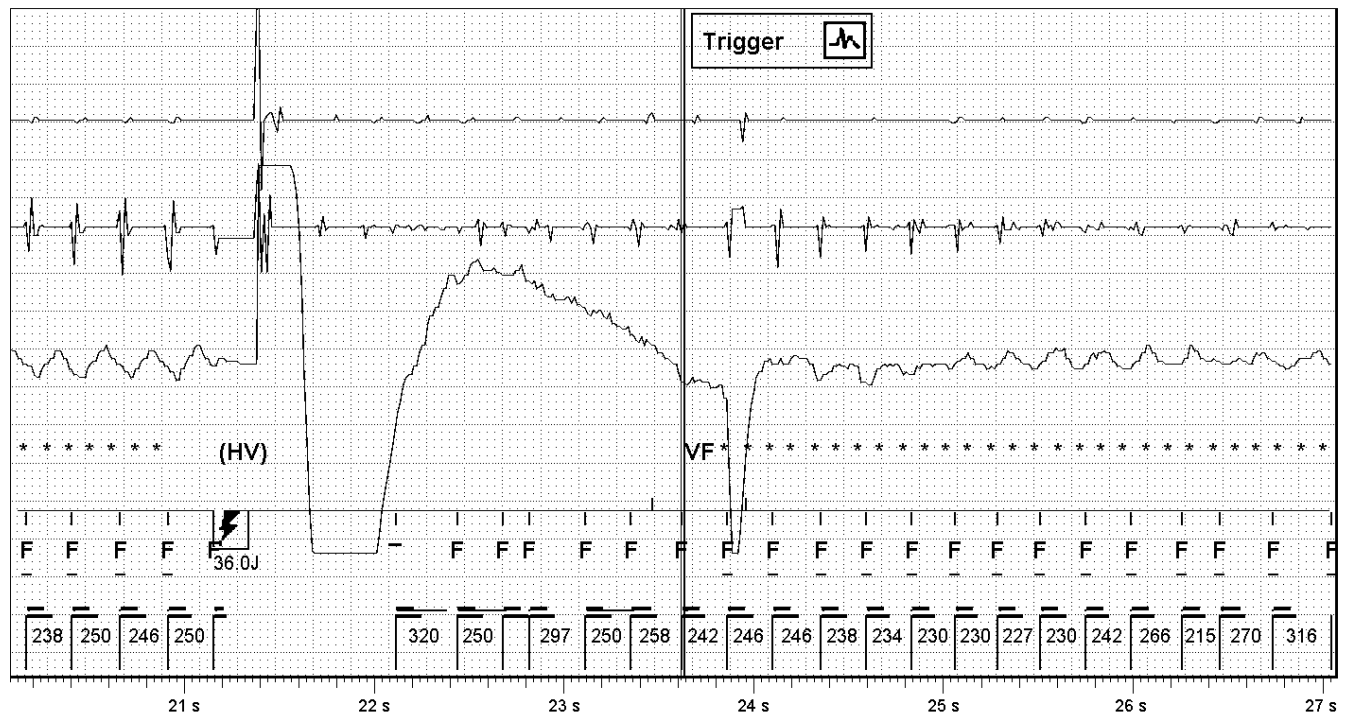


Episode: VF (240 bpm / 250 ms) (Continued)

VT/VF Episode 6 of 6

Page 4 of 5

9 Mar 2022 13:21

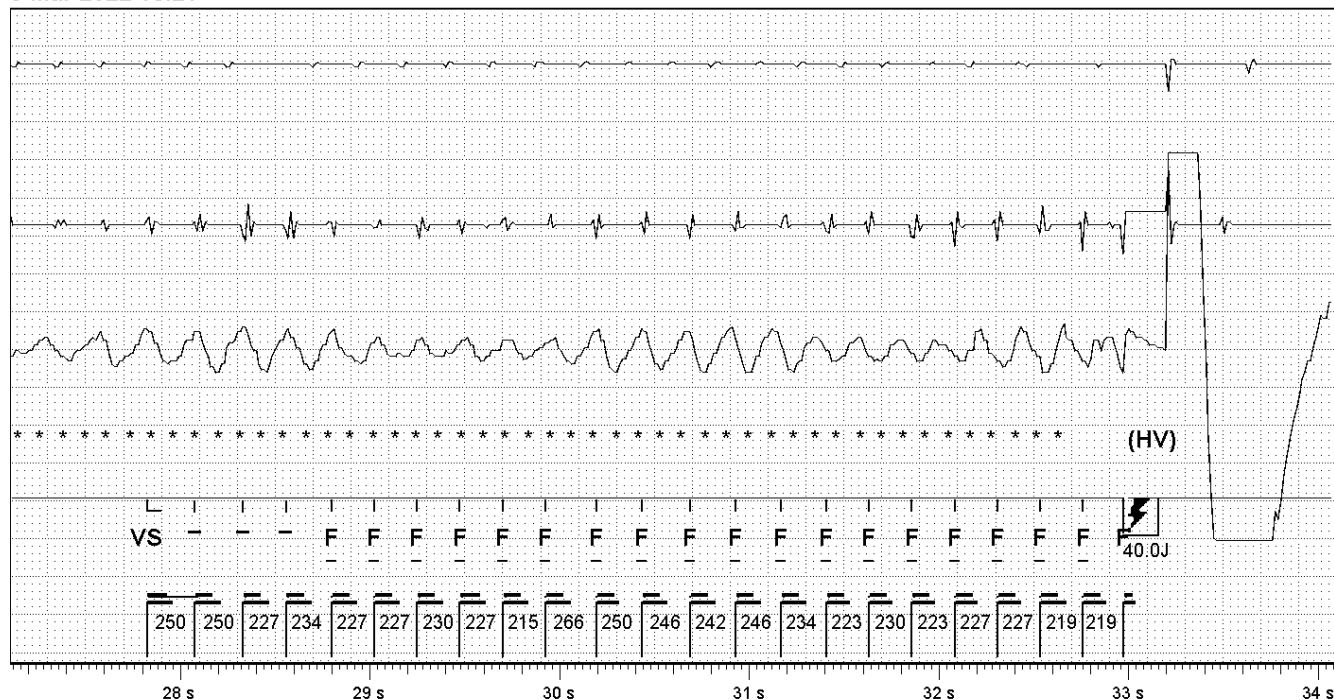


1: A Sense Amp AutoGain (3,4 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,7 mm/mV)

3: Discrimination AutoGain (2,6 mm/mV)

Sweep Speed: 25 mm/s

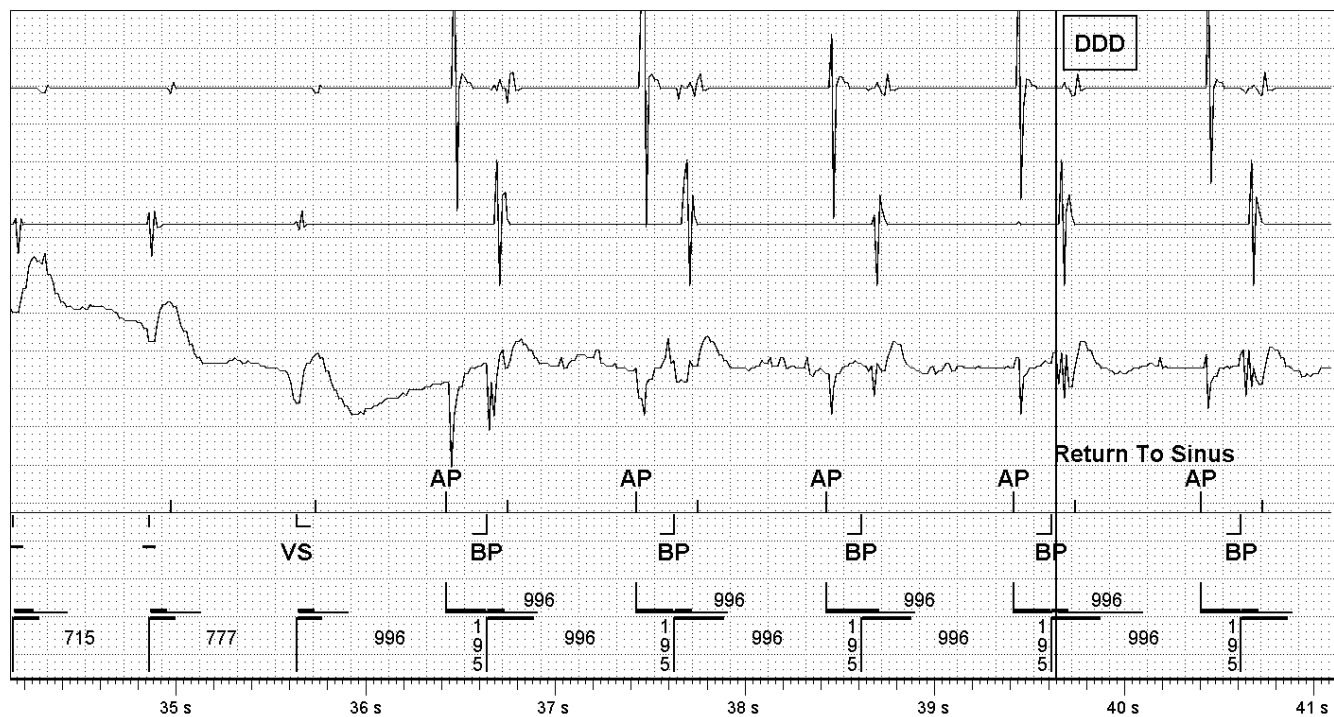


Episode: VF (240 bpm / 250 ms) (Continued)

VT/VF Episode 6 of 6

Page 5 of 5

9 Mar 2022 13:21

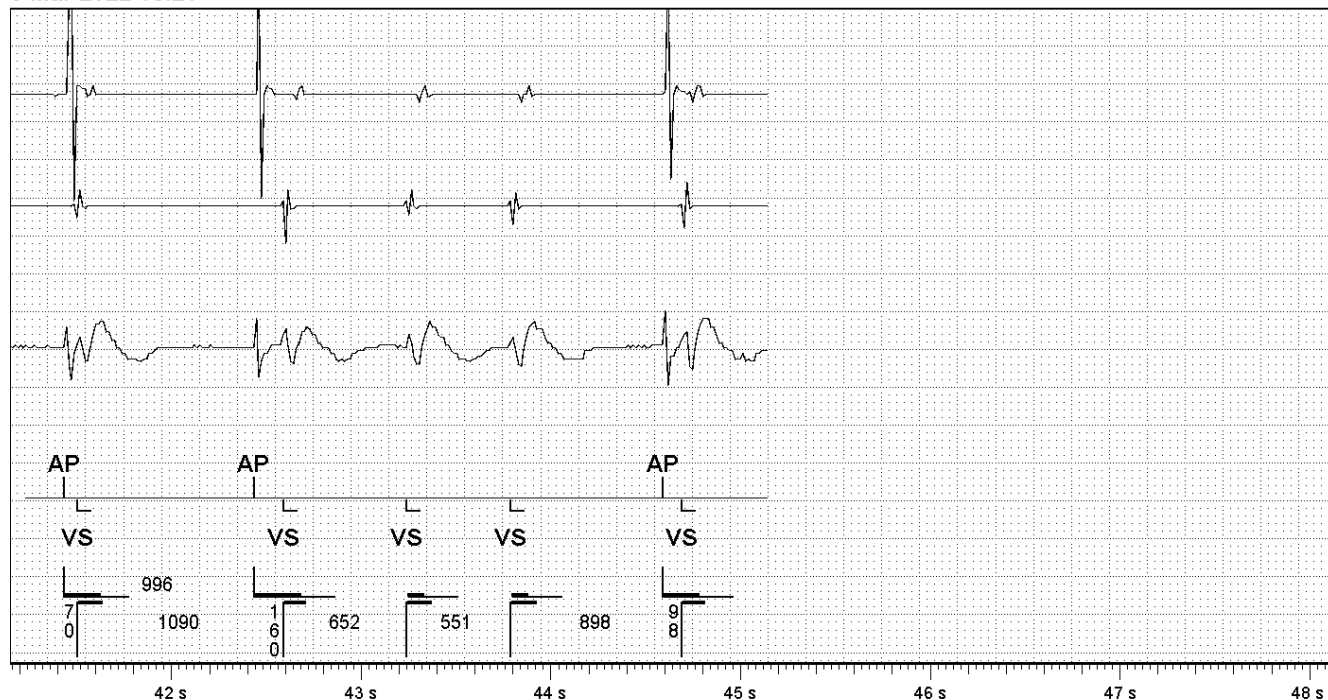

1: A Sense Amp AutoGain (3,4 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,7 mm/mV)

3: Discrimination AutoGain (2,6 mm/mV)

Sweep Speed: 25 mm/s

**Episode: VF (235 bpm / 255 ms)**

VT/VF Episode 5 of 6

Page 1 of 4

9 Mar 2022 13:21

Duration 00:17 (M:S)

Detection Criteria 214 bpm

Alerts



Three or more VT/VF episodes in 24 hours

Therapy1) ATP x 1
2) Defib 36,0 J (844V)**Results**1) VF
2) Below Rate Detection
(CL 860 ms)**HV Therapy**Last HV Lead Impedance 40 Ω
Last Charge Time 6,6 sec
Delivered PW +5,0 ms, -5,0 ms**ATP Therapy Details****VF ATP Therapy****Therapy**

Burst 1 216, 216, 216, 216, 216, 216 ms

Diagnosis Summary

Diagnosis

VF

Time to Diagnosis

3,25 sec

Rate (CL)

235 bpm (255 ms)

DeFT Response™ Settings (Shock Waveform)

Waveform Biphasic

Waveform Mode Tilt

Shock Configuration RV to SVC & Can

RV Polarity Anode(+)

Defib Max Sensitivity 0.5 mV

VF Shocks (Defib)

VT Shocks (CVRT)

1st Phase

65 %

Same as Defib

2nd Phase

65 %


Same as Defib

Episode: VF (235 bpm / 255 ms)

VT/VF Episode 5 of 6

Page 2 of 4

9 Mar 2022 13:21

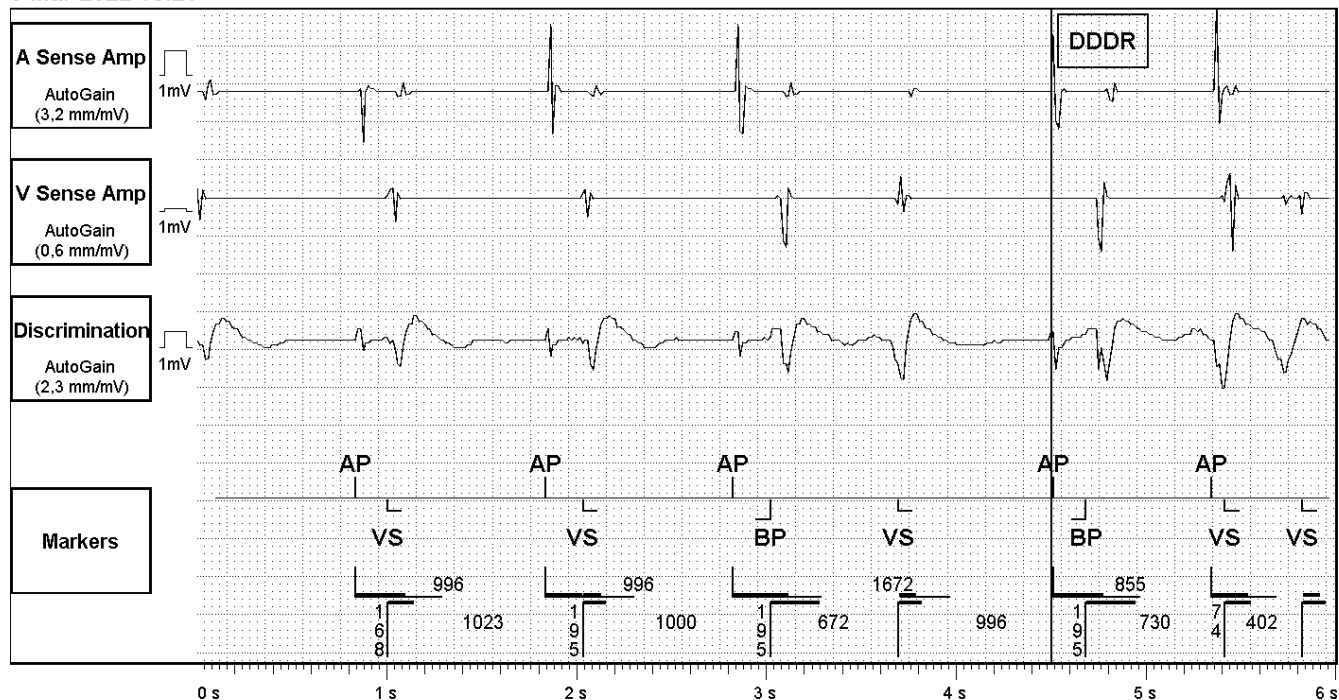


1: A Sense Amp AutoGain (3,2 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,6 mm/mV)


3: Discrimination AutoGain (2,3 mm/mV)

Sweep Speed: 25 mm/s

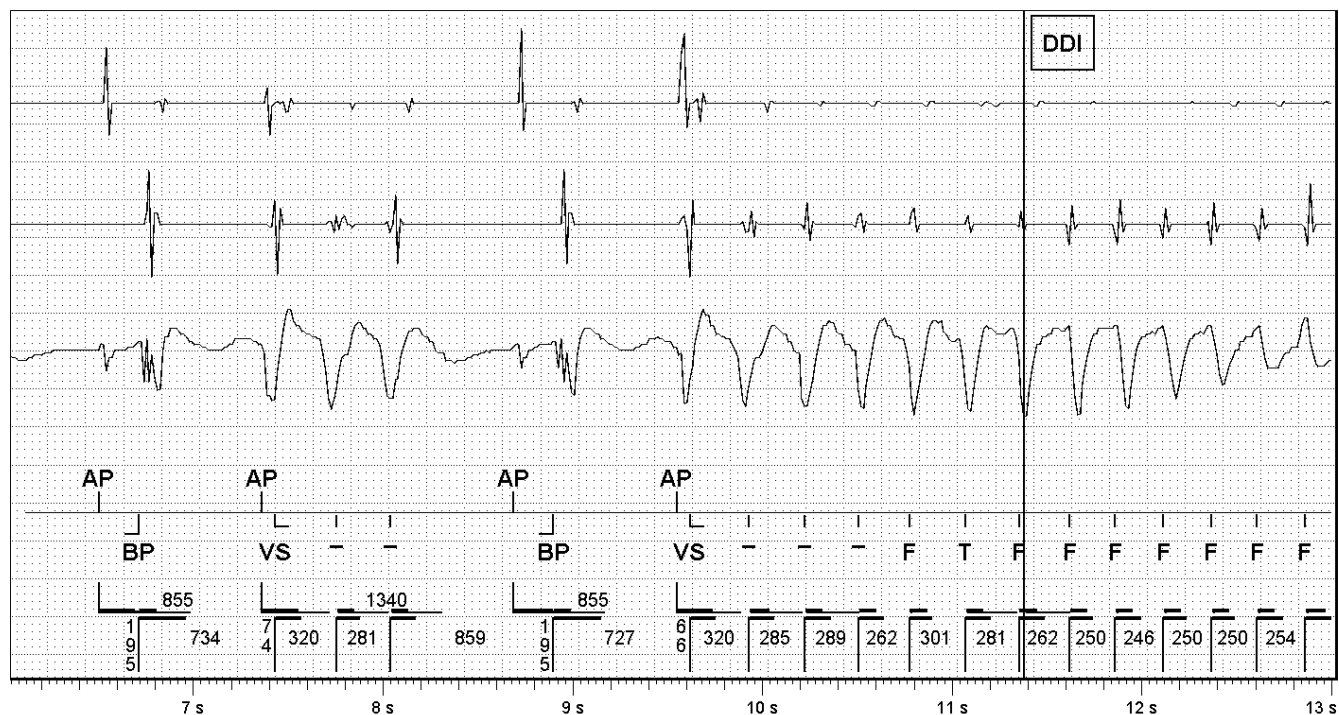

Page 3 of 4

Sweep Speed: 25 mm/s

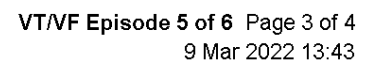


Episode: VF (235 bpm / 255 ms) (Continued)

VT/VF Episode 5 of 6

Page 4 of 4

9 Mar 2022 13:21

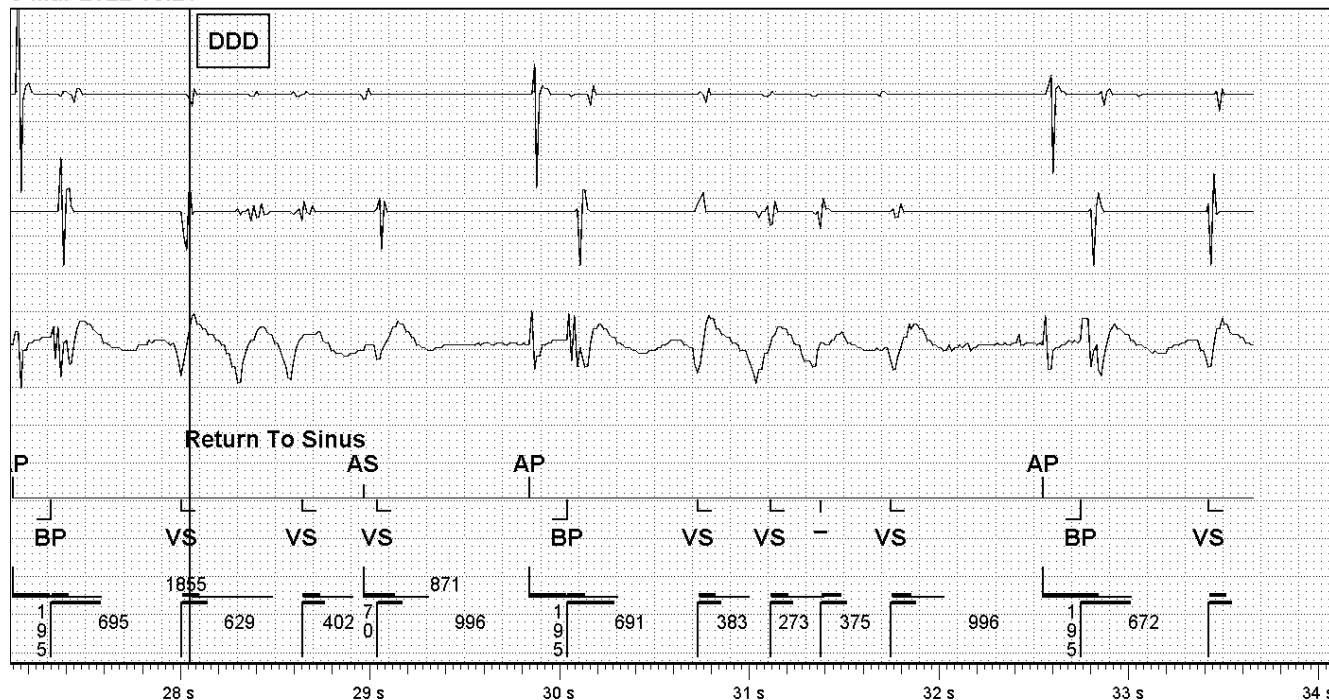

1: A Sense Amp AutoGain (3,2 mm/mV)
2: V Sense Amp AutoGain (0,6 mm/mV)
3: Discrimination AutoGain (2,3 mm/mV)

4: Markers

Sweep Speed: 25 mm/s

**Episode: VF (235 bpm / 255 ms)**

VT/VF Episode 4 of 6

Page 1 of 4

9 Mar 2022 13:20

Duration 00:17 (M:S)
Detection Criteria 214 bpm

Alerts



Three or more VT/VF episodes in 24 hours

Therapy1) ATP x 1
2) Defib 36,0 J (844V)**Results**1) VF
2) Below Rate Detection
(CL 880 ms)**HV Therapy**Last HV Lead Impedance 40 Ω
Last Charge Time 7,5 sec
Delivered PW +4,8 ms, -4,8 ms**ATP Therapy Details****VF ATP Therapy****Therapy**

Burst 1 216, 216, 216, 216, 216, 216 ms

Diagnosis Summary

Diagnosis


VFTime to Diagnosis 2,75 sec
Rate (CL) 235 bpm (255 ms)**DeFT Response™ Settings (Shock Waveform)**Waveform Biphasic
Waveform Mode Tilt
Shock Configuration RV to SVC & Can
RV Polarity Anode(+)
Defib Max Sensitivity 0.5 mVVF Shocks (Defib)
VT Shocks (CVRT)**1st Phase**65 %
Same as Defib**2nd Phase**65 %
Same as Defib

Episode: VF (235 bpm / 255 ms)

VT/VF Episode 4 of 6

Page 2 of 4

9 Mar 2022 13:20

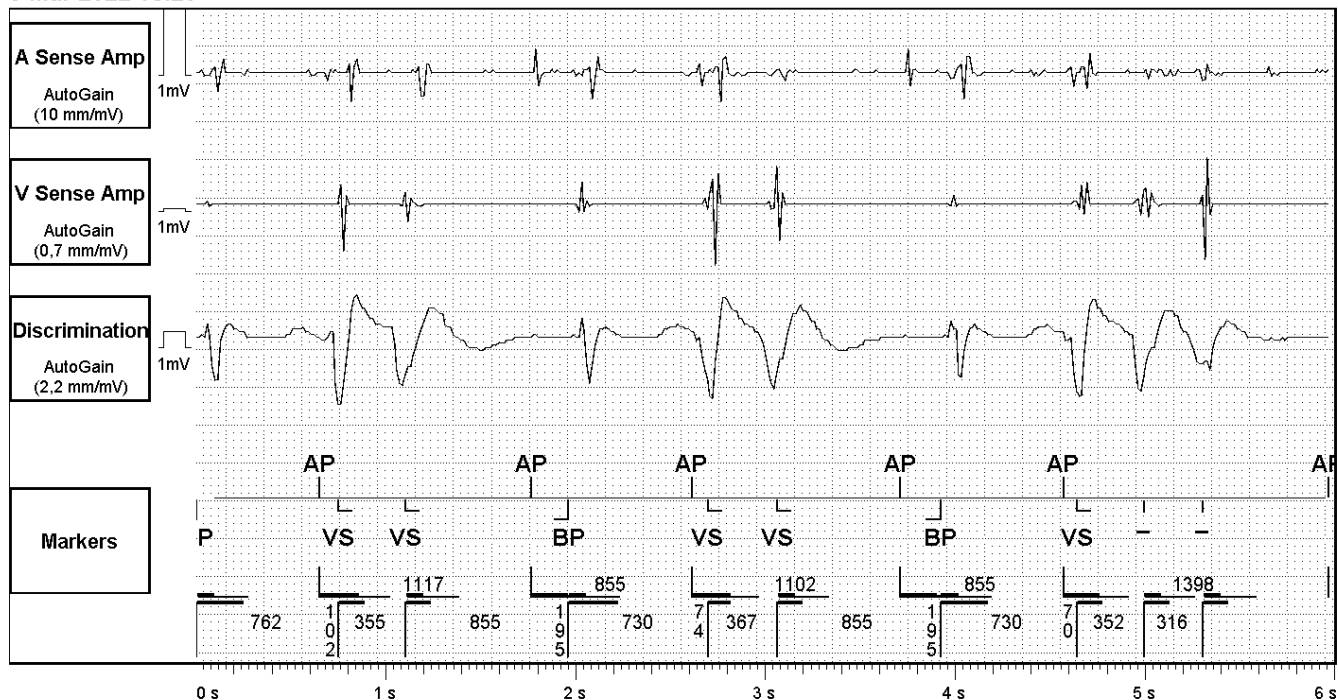


1: A Sense Amp AutoGain (10 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0.7 mm/mV)

3: Discrimination AutoGain (2.2 mm/mV)

Sweep Speed: 25 mm/s

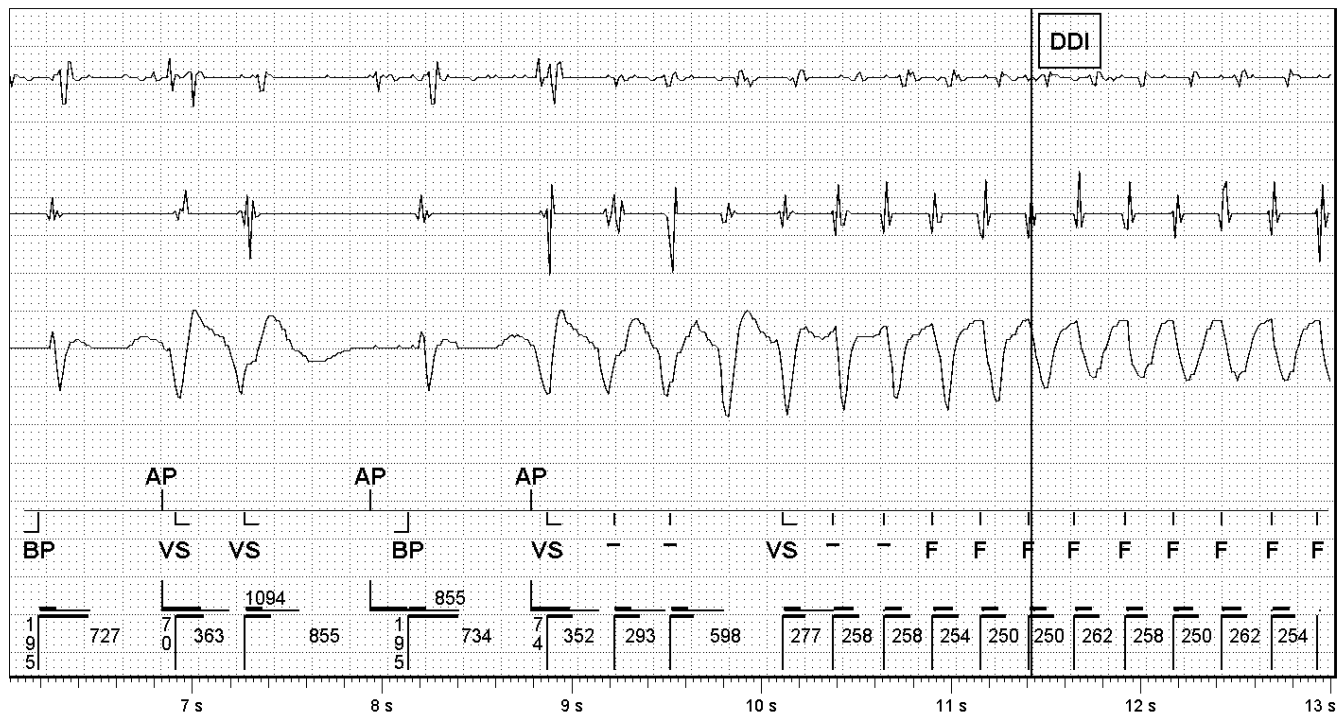


Episode: VF (235 bpm / 255 ms) (Continued)

VT/VF Episode 4 of 6

Page 3 of 4

9 Mar 2022 13:20

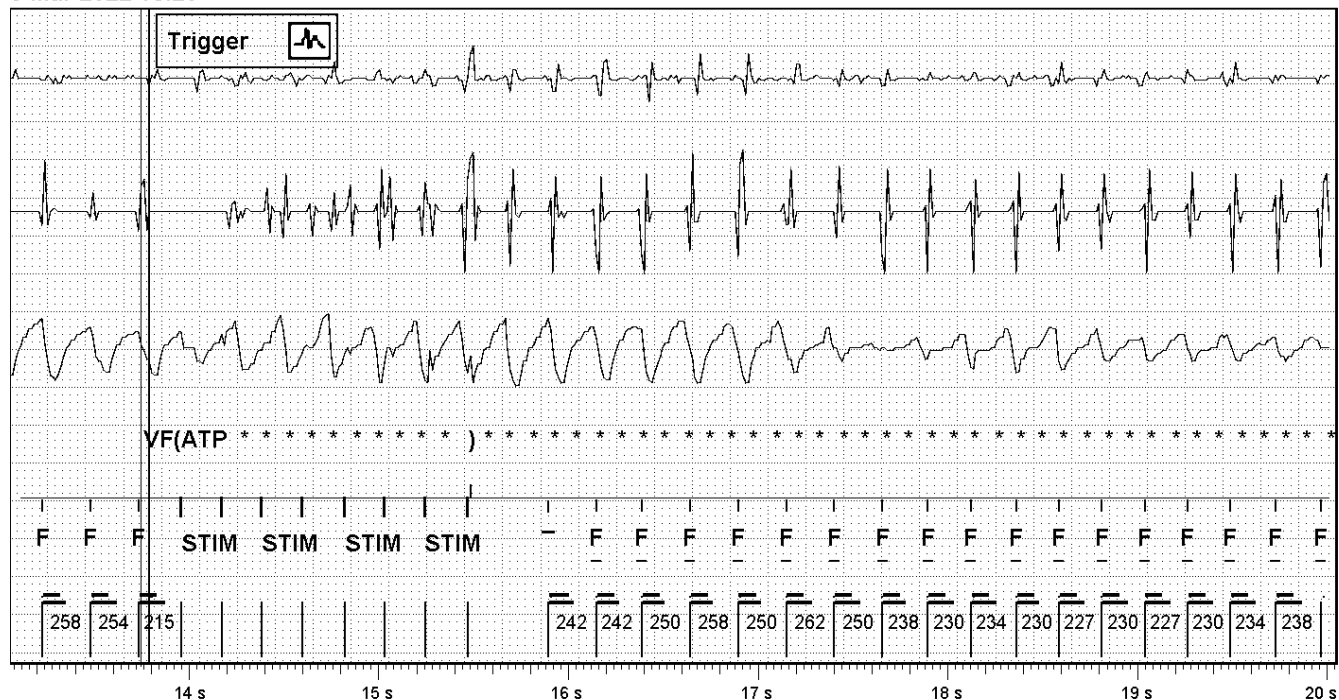


1: A Sense Amp AutoGain (10 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,7 mm/mV)

3: Discrimination AutoGain (2,2 mm/mV)

Sweep Speed: 25 mm/s

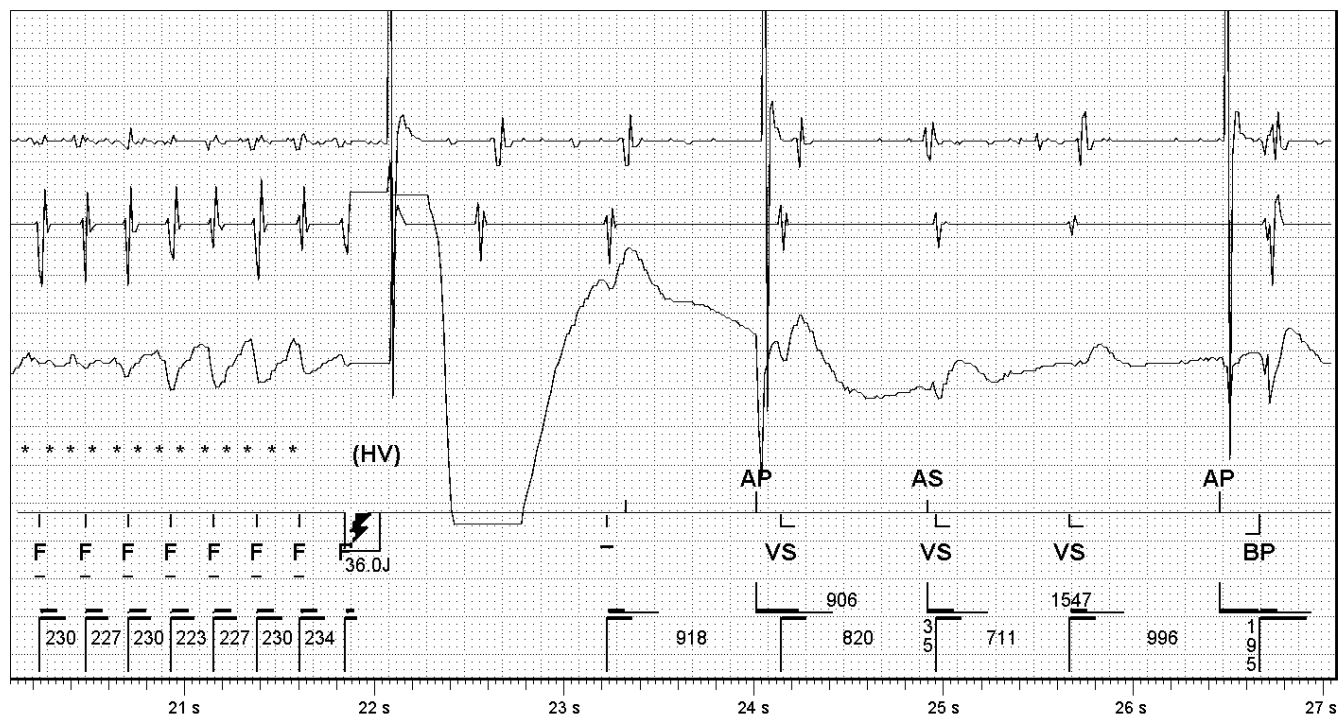


Episode: VF (235 bpm / 255 ms) (Continued)

VT/VF Episode 4 of 6

Page 4 of 4

9 Mar 2022 13:20

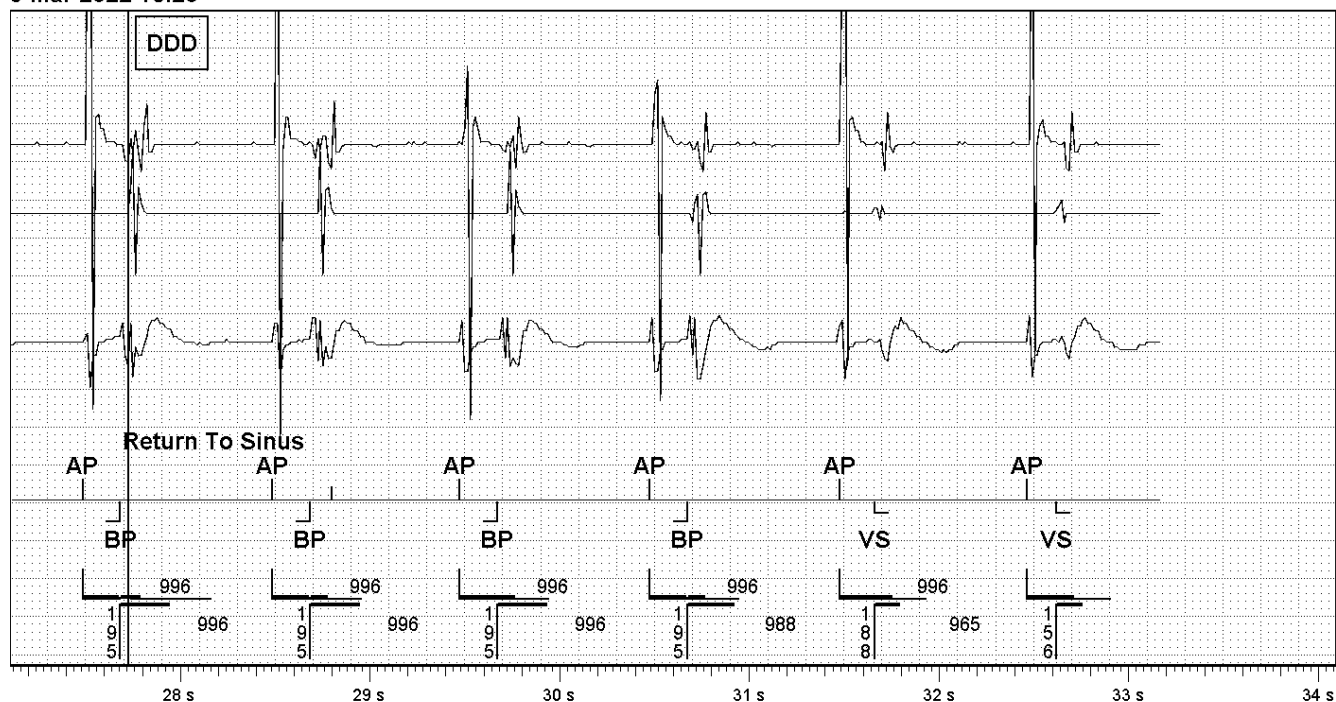

1: A Sense Amp AutoGain (10 mm/mV)

4: Markers

2: V Sense Amp AutoGain (0,7 mm/mV)

3: Discrimination AutoGain (2,2 mm/mV)

Sweep Speed: 25 mm/s

Test Results: Atrium

Page 1 of 4

Atrial Capture Test

No result recorded
0,75 V @ 0,7 ms (Bi)

9 Mar 2022
9 Mar 2022


Atrial Sense Test

No result recorded
None (Bi)

9 Mar 2022
9 Mar 2022


ACap™ Confirm Trend

Bi

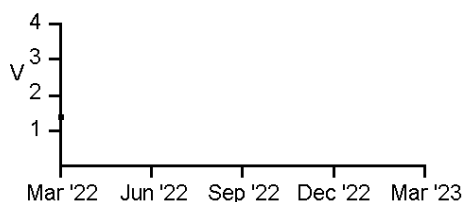

Atrial Sense Amplitude Trend


Bi

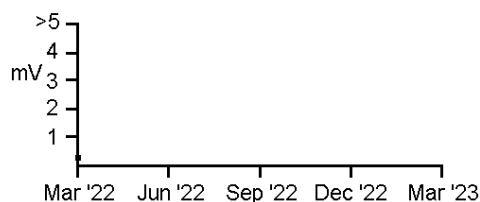

Atrial Lead Impedance:


No recent measurement

1-year trend

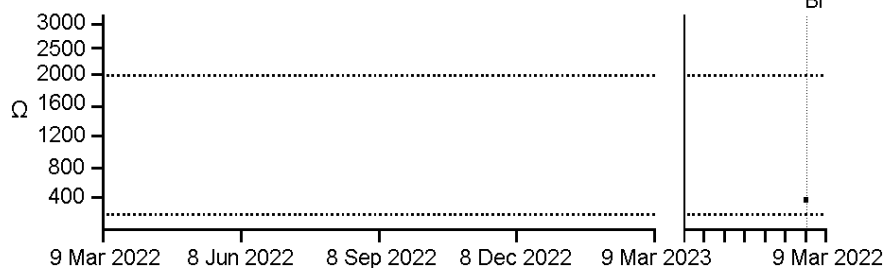

Last 7 Days

Bi

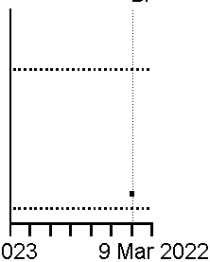

Programmed Pulse Configuration	Bipolar
Lead Monitoring	Monitor
Lower Limit	200 Ω
Upper Limit	2.000 Ω

Configurations

Bipolar

First Measurement

440 Ω (8 Mar 2022)

Lifetime Range

390 - 440 Ω



Test Results: Right Ventricle

Page 2 of 4

Right Ventricular Capture Test

No result recorded

0,5 V @ 0,5 ms (Bi)

9 Mar 2022

9 Mar 2022

Ventricular Sense Test

No result recorded


5,6 - 6,0 mV (RV Bi)

9 Mar 2022

9 Mar 2022

RV Sense Amplitude Trend


RV Bi

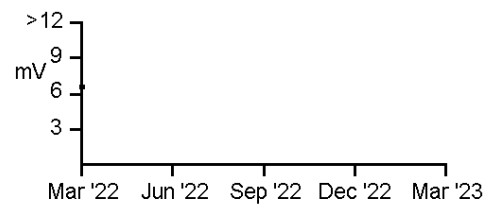

Right Ventricular Lead Impedance:


No recent measurement

1-year trend

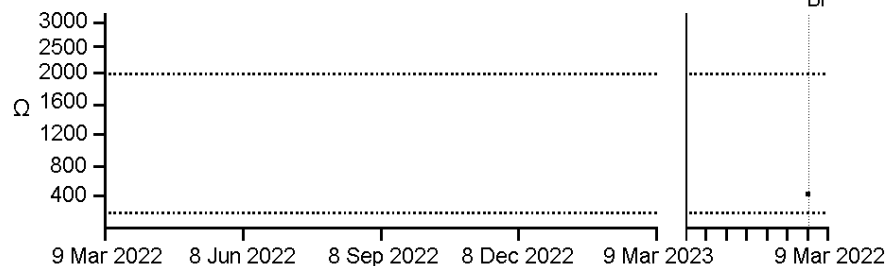

Last 7 Days

Bi

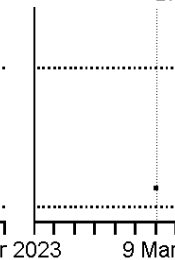

Programmed Pulse

Configuration

Lead Monitoring

Lower Limit

Upper Limit

Bipolar

Monitor

200 Ω

2.000 Ω

Configurations

Bipolar

First Measurement

510 Ω (8 Mar 2022)

Lifetime Range

410 - 510 Ω

Test Results: Left Ventricle

Page 3 of 4

Left Ventricular Capture Test

No result recorded
1,0 V @ 0,5 ms (D1 - M2)

9 Mar 2022
9 Mar 2022

LV2 Capture Test

No result recorded
No previous results

9 Mar 2022

LV Lead Impedance:

No recent measurement


LV2 Lead Impedance:

380 Ω (Mid 3 - RV Coil)

9 Mar 2022

LV: 1-year trend

Last 7 Days

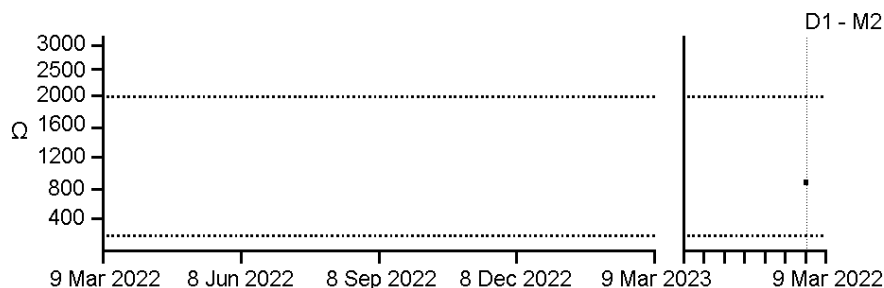


LV Pulse Configuration	D1 - M2
LV2 Pulse Configuration	M3 - RVc
Lead Monitoring	Monitor
Lower Limit	200 Ω
Upper Limit	2.000 Ω

Configurations	First Measurement	Lifetime Range	Configurations	First Measurement	Lifetime Range
Distal tip 1 - Mid 2	1150 Ω (8 Mar 2022)	750 - 1150 Ω	Mid 3 - RV Coil	430 Ω (8 Mar 2022)	380 - 450 Ω

Test Results: HV Lead Impedance and Battery

Page 4 of 4

HV Lead Impedance:

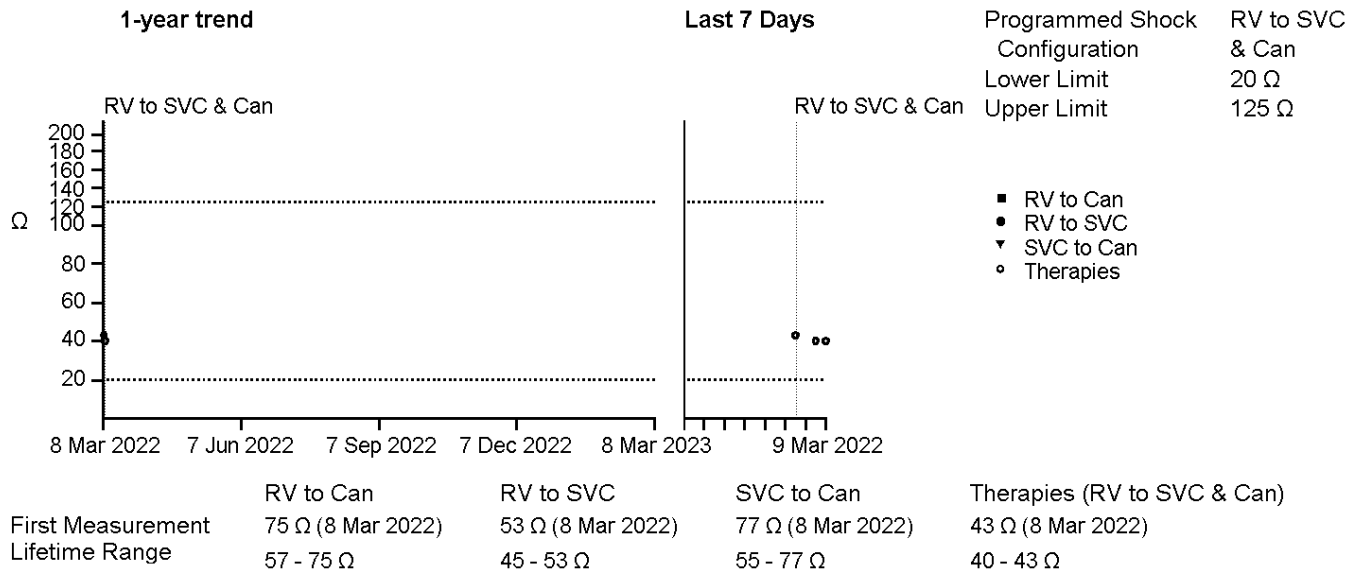 40 Ω (RV to SVC & Can) 9 Mar 2022


Battery Information

Longevity: 5, 1-5, 4 yrs



Last Max Charge

8,6 sec (Therapy)

Longevity estimate is based on patient history

Battery Current

 21 μ A

Remaining Capacity to ERI

91%

Presenting Rhythm Freeze

Page 1 of 1

9 Mar 2022 13:37

Key Parameters

Mode	DDDR
Base Rate	70 bpm
Rest Rate	60 bpm
Paced AV Delay	200 ms
Sensed AV Delay	150 ms
Max Track Rate	110 bpm
Max Sensor Rate	110 bpm
Hysteresis Rate	Off
ACap™ Confirm	Monitor
RVCap™ Confirm	Off
LVCap™ Confirm	Off
AF Suppression™	Off
SyncAV™ CRT Delta	Off
Rate Responsive AV Delay	Medium
Rate Responsive PVARP/V Ref	Low
Ventricular Safety Standby	On

1: II AutoGain (100,0 mm/mV)


2: Markers

3: A Sense Amp AutoGain (10,0 mm/mV)

4: RV Bipolar AutoGain (4,4 mm/mV)

5: LV Distal tip 1 - Mid 2 AutoGain (1,1 mm/mV)

Sweep Speed: 25 mm/s

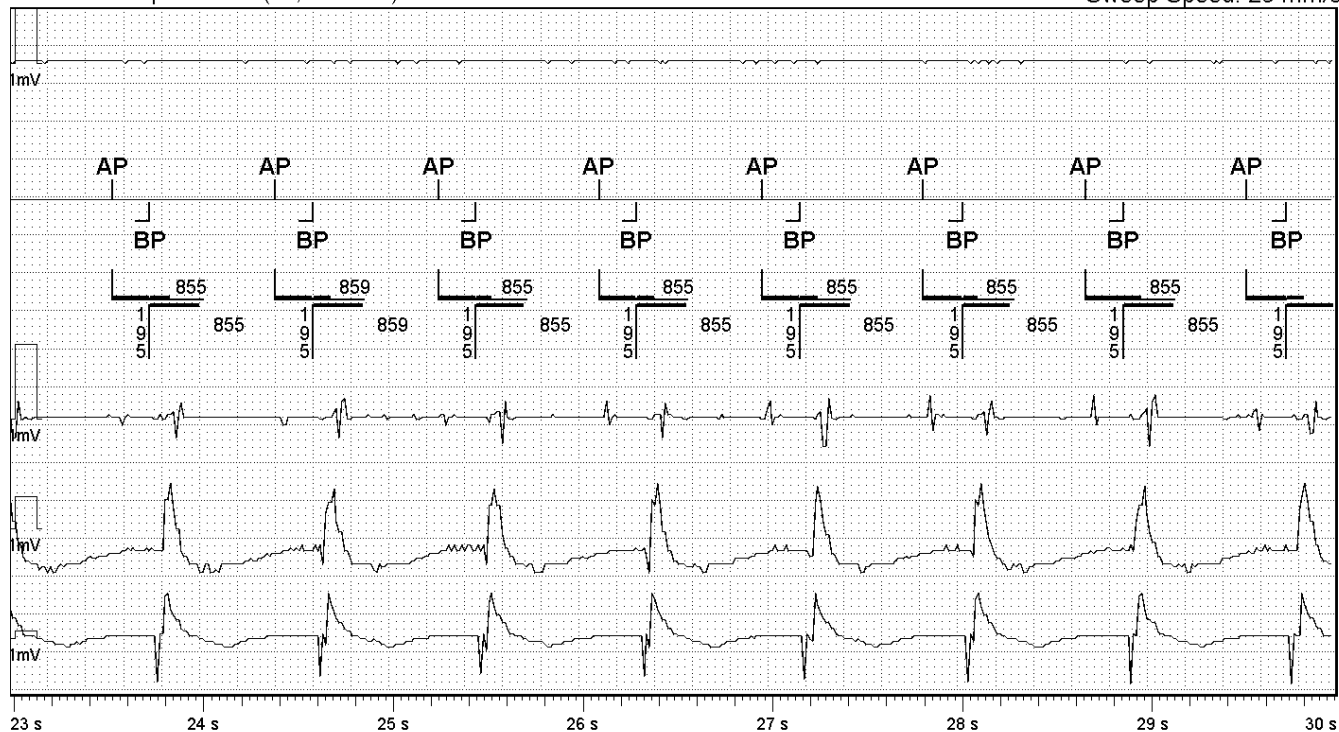

Wrap-up™ Overview

 **2 Alerts**

Page 1 of 1

Patient


Date of Birth 8 Jun 1944
EF % Unknown

Indications for Implant

Device	Manufacturer	Model	Serial	Implant Date
CRT-D	St. Jude Medical	Quadra Assura MP™ 3371-40QC	9973677	8 Mar 2022
A Lead				
RV Lead	St. Jude Medical	Durata™ 7120Q / 58 cm	BNP078743	8 Mar 2022
LV Lead	St. Jude Medical	Quartet™ 1457Q / 86 cm	DBU030198	8 Mar 2022

Battery


Longevity: 5,1-5,4 yrs

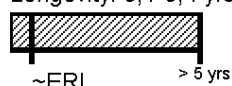
Last Max Charge 8,6 sec (Therapy)
Battery Current 21 uA
Remaining Capacity to ERI 91%

Longevity estimate is based on patient history

Test Results 9 Mar 2022

 Automatic

	Capture	Sense	Lead Impedance
A	Not Performed	Not Performed	Not Performed
RV	Not Performed	Not Performed	Not Performed
LV	Not Performed		Not Performed
HV			40 Ω (RV to SVC & Can) 

Parameters

Mode	DDDR	Zone Configuration	VT	VF
Base Rate	▶90 bpm	Detection Criteria	171 bpm	214 bpm
Max Track Rate	110 bpm	Therapy (ENABLED)	ATP x3	ATP x1
Paced AV Delay	200 ms		25,0 J	36,0 J
Sensed AV Delay	150 ms		36,0 J	40,0 J
Ventricular Pacing	LV → RV, 30 ms		40,0 J x2	40,0 J x4

Programming Changes

	Initial	Final
AMS Base Rate	70 bpm	↳90 bpm
Base Rate	70 bpm	▶90 bpm
Rest Rate	60 bpm	▶Off

Alerts

Episodes with alert conditions (3)
A. Pulse Amplitude Safety Margin < 2:1

- ▶ Manual-programmed
- ↳ Auto-programmed