

ZOOM ® View™ Report Created 04 Jun 2018

Combined Follow-up Report Hucík, Stepán

Date of Birth 7 Jan 1944

INCEPTA CRT-D P162/119874 Device

Tachy Mode Monitor + Therapy Last Office Interrogation

19 Jun 2017 Implant Date 12 Mar 2013

My Alerts

a 04 Jun 2018 09:53 Check Atrial and LV Leads

Events Since Last Reset (14 Mar 2013)

VT at 200 min⁻¹, ATPx4, 17J 30 Aug 2014 13:14 VT at 181 min⁻¹, **♦** ATPx4, 17J 09 Aug 2014 16:46

See last page for full list of Events Since Last Reset.

Battery OK

> Approximate time to explant: 3.5 years Charge Time 10.7 s

> Last Capacitor Re-form 27 Apr 2018 11:27

One Year Remaining **Explant**

Leads Data Most Recent Implant Previous Session 12 Mar 2013 Atrial N/R mV N/R mV Intrinsic Amplitude N/R mV Pace Impedance $N/R \Omega$ >2500 Ω >2500 Ω N/R V @ N/R ms Pace Threshold 4.0 V @ 0.4 ms **Right Ventricular** Intrinsic Amplitude 7.0 mV PVC mV 6.0 mV @ 44 min⁻¹ Pace Impedance 650Ω 550Ω 507 Ω Pace Threshold 0.5 V @ 0.4 ms 0.7 V @ 0.4 ms 1.1 V @ 0.4 ms Left Ventricular PVC mV Intrinsic Amplitude 12.0 mV 20.7 mV Pace Impedance 650Ω 544Ω 545 Ω Pace Threshold 1.3 V @ 0.4 ms 1.1 V @ 0.4 ms 1.2 V @ 0.4 ms **Shock Vector** Shock Impedance 45 Ω 56 Ω 52 Ω

Settings

Ventricular Tachy Settings

31J, 41J, 41Jx6 VF 230 min⁻¹ ATP 205 min⁻¹ Scan VT 31J, 41J, 41Jx4 Ramp 41J, 41J, 41Jx3 VT-1 160 min⁻¹ Scan Ramp

Brady/CRT Settings

Mode **VVIR** Pacing Output Right Ventricular Lower Rate Limit 60 min⁻¹ 3.0 V @ 0.4 ms Maximum Sensor Rate 130 min⁻¹ Left Ventricular 3.0 V @ 0.4 ms RV-Refractory (RVRP) 230 - 250 ms Sensitivity LV-Refractory (LVRP) 250 ms Atrial AGC 0.25 mV Ventricular Pacing Chamber BiV Right Ventricular AGC 0.6 mV LV Offset AGC 1.0 mV Left Ventricular 0 ms

Leads Configuration (Pace/Sense)

Bipolar/Off Atrial Right Ventricular **Bipolar** Left Ventricular LVring>>RV / LVtip>>LVring


Sensors


Accelerometer On

2868 Software Version: 4.04 P162 Firmware Version: B_v1.02.00(4.01)

© 2016 **Boston Scientific Corporation** or its affiliates. All rights reserved. Page 1 of 6


Clinician Signature:

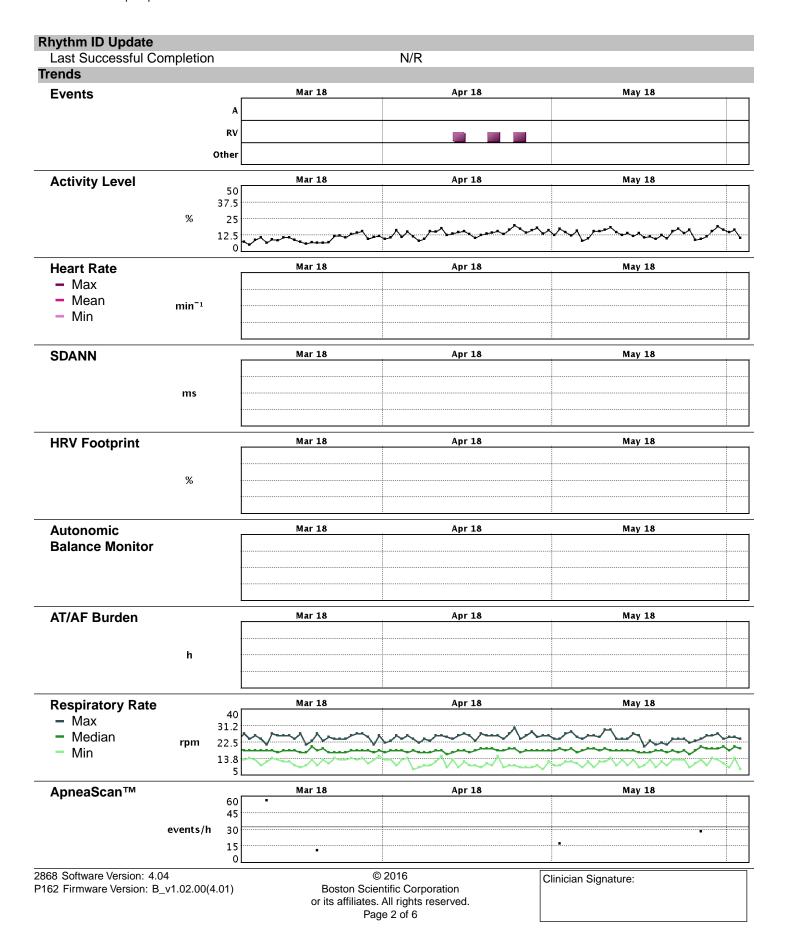


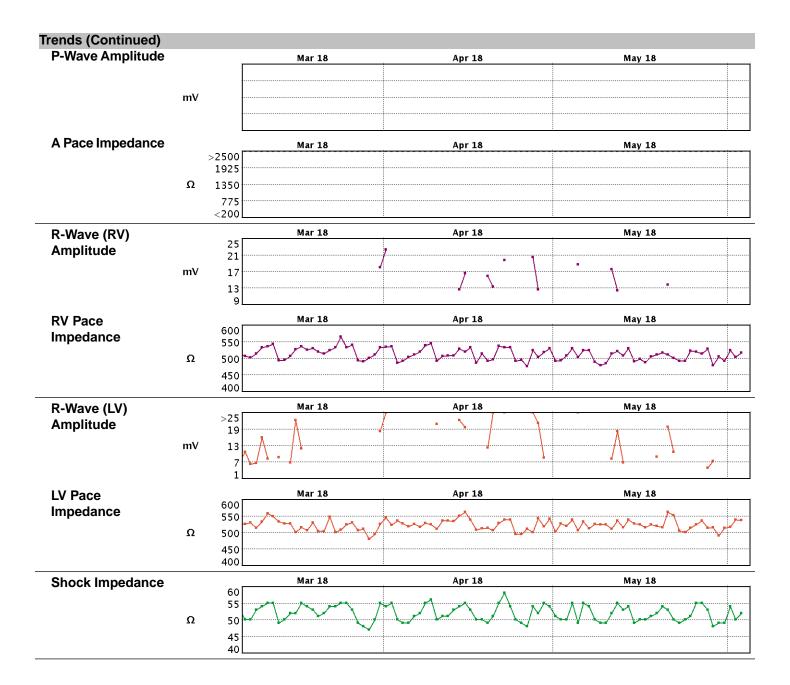


Ventricular Tachy Counters	Since Last Reset 14 Mar 2013	Device Totals
Ventricular Episode Counters		
Total Episodes	7.2K	7.2K
Treated		
VF Therapy	0	0
VT Therapy	4	4
VT-1 Therapy	17	17
Commanded Therapy	0	0
Non-Treated		
No Therapy Programmed	0	0
Nonsustained Episodes	7.2K	7.2K
Other Untreated Episodes	7	7
Ventricular Therapy Counters		
ATP Delivered	29	29
ATP % Successful	67	67
Shocks Delivered	6	6
First Shock % Successful	100	100
Shocks Diverted	2	2
Brady/CRT Counters	Reset Before Last	Since Last Reset
	2 day(s)	1907 day(s)
	12 Mar 2013 to 14 Mar 2013	3 14 Mar 2013 to Today
Counters		
% A Paced	100	0
% RV Paced	90	90
% LV Paced	99	97
Intrinsic Promotion		
Rate Hysteresis		
% Successful	0	0
Atrial Arrhythmia		
% AT/AF	N/R	0
(One Year of Data)		
Total Time in AT/AF (min)	N/R	0.0
(One Year of Data)		
Episodes by Duration		
< 1 minute	0	0
1 min - < 1 hr	0	0
1 hr - < 24 hr	0	0
24 hr - < 48 hr	0	0
> 48 hr	0	0
Total PACs	0	0
Ventricular Counters		
Total PVCs	14	88
Three or More PVCs	1	8


Histograms Reset Before Last (12 Mar 2013)

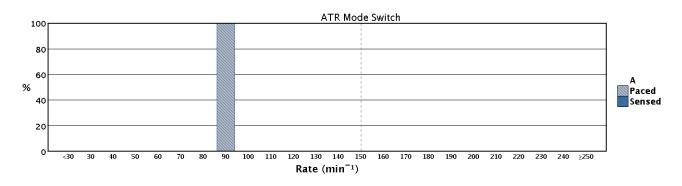

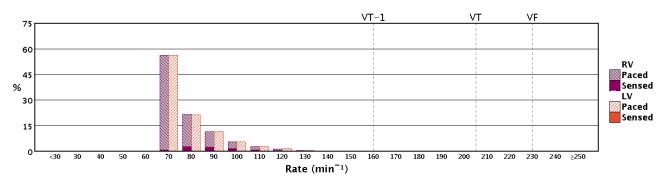

Histograms Since Last Reset (14 Mar 2013)

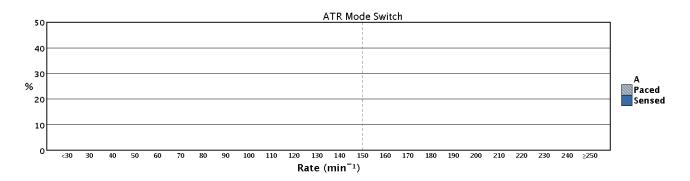

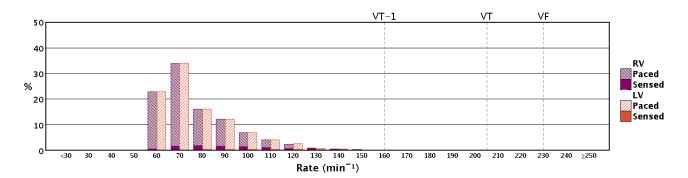

2868 Software Version: 4.04 P162 Firmware Version: B_v1.02.00(4.01) © 2016

Boston Scientific Corporation
or its affiliates. All rights reserved.

Page 5 of 6

Clinician Signature:

Heart Rate Variability

Last Measured:

Last Measured Data unavailable

Reference:

Reference Data Unavailable

All Events Since Last Reset (14 Mar 2013)

I Events Since Last Reset	(14 Mar 2013)
25 Apr 2018 08:08	NonSustV at 169 min ⁻¹ , Nonsustained
20 Apr 2018 18:38	NonSustV at 153 min ⁻¹ , Nonsustained
14 Apr 2018 18:31	NonSustV at 156 min ⁻¹ , Nonsustained
14 Apr 2018 17:08	NonSustV at 159 min ⁻¹ , Nonsustained
14 Apr 2018 17:02	NonSustV at 159 min ⁻¹ , Nonsustained
14 Apr 2018 17:01	NonSustV at 158 min ⁻¹ , Nonsustained
14 Apr 2018 17:00	NonSustV at 161 min ⁻¹ , Nonsustained
14 Apr 2018 16:55	NonSustV at 164 min ⁻¹ , Nonsustained
14 Apr 2018 16:54	NonSustV at 157 min ⁻¹ , Nonsustained
14 Apr 2018 14:00	NonSustV at 161 min ⁻¹ , Nonsustained
03 Dec 2015 12:23	VT-1 at 181 min ⁻¹ , No Therapy
03 Dec 2015 12:11	VT-1 at 172 min ⁻¹ , No Therapy
03 Dec 2015 11:32	VT-1 at 146 min ⁻¹ , No Therapy
19 Oct 2015 14:51	VT-1 at 175 min ⁻¹ , No Therapy
23 Sep 2015 16:09	VT-1 at 163 min ⁻¹ , ATPx1 VT-1 at 173 min ⁻¹ , ATPx1
23 Sep 2015 15:54 23 Sep 2015 14:11	VT-1 at 173 min 7, ATPX1 VT-1 at 157 min -1, ATPX1
•	VT-1 at 162 min ⁻¹ , ♦ ATPx4
23 Sep 2015 13:11	
23 Sep 2015 13:10	VT-1 at 179 min ⁻¹ , ♦ ATPx5
23 Sep 2015 13:07	VT-1 at 173 min ⁻¹ , ♦ ATPx2
21 Sep 2015 20:26	VT-1 at 178 min ⁻¹ , ATPx1
21 Sep 2015 19:19	VT-1 at 186 min ⁻¹ , ATPx1
17 Sep 2015 12:20	VT-1 at 181 min ⁻¹ , ♦ ATPx8, 41J
17 Sep 2015 12:18	VT-1 at 181 min ⁻¹ , ♦ ATPx5
17 Sep 2015 10:44	VT-1 at 196 min ⁻¹ , ATPx1
30 Aug 2015 17:02	VT-1 at 155 min⁻¹, � ATPx3
30 Aug 2015 16:19	VT-1 at 165 min ⁻¹ , No Therapy
29 Aug 2015 14:33	VT-1 at 170 min ⁻¹ , No Therapy
29 Aug 2015 14:32	VT-1 at 186 min ⁻¹ , ATPx1
29 Aug 2015 14:20	VT-1 at 176 min ⁻¹ , ATPx1
18 May 2015 10:09	VT-1 at 180 min⁻¹, No Therapy
25 Oct 2014 11:49	VT-1 at 183 min ⁻¹ , ♦ ATPx8, 41J
25 Oct 2014 11:38	VT-1 at 155 min ⁻¹ , ATPx1
10 Oct 2014 12:13	VT-1 at 168 min ⁻¹ , ♦ ATPx8, 41J
30 Aug 2014 13:14	VT at 200 min⁻¹, � ATPx4, 17J
09 Aug 2014 16:46	VT at 181 min⁻¹, � ATPx4, 17J
16 May 2013 17:54	VT at 178 min⁻¹, � ATPx2
17 Apr 2013 15:14	VT at 179 min ⁻¹ , ATPx4, 17J
•	, - , -

2868 Software Version: 4.04 P162 Firmware Version: B_v1.02.00(4.01) © 2016

Boston Scientific Corporation
or its affiliates. All rights reserved.

Page 6 of 6

Clinician Signature: