Scientific Scientific

ZOOM ® View™		Report Created 11 Dec 2020	
Combined Follow-up Report		Last Office Interrogation	
		10-JUL-2013	
		Implant Date	
Device	ALTRUA 50 S508/576134	N.R. N.R. N.R.	

Events Since Last Reset

09-JUL-2013 12:21 Tach at 190 min⁻¹
09-JUL-2013 12:14 Tach at 280 min⁻¹
See last page for full list of Events Since Last Reset.

Battery	1	Go	000	l

Longevity Remaining 3.5 years

@ Current Pacing Percentage(s)

Magnet Rate 100 min⁻¹

Elective Replacement (ERT)	/	\
End of Life (EOL)		
		_

Beginning of Life (BOL)

	Lind of Line (LOL)			
Leads Data	Implant	Previous Session	Present Session	
	N.R. N.R. N.R.		10-JUL-2013	
Intrinsic Amplitude	N.R. mV	5 mV	3.3 mV	
Pace Impedance	N.R. Ω	$580~\Omega$	$560~\Omega$	
Pace Threshold	N.R. V		1.4 V @ 0.6 ms	

Settings


Tachy Settings

Tachy Detection On
Rate 170 min⁻¹
Duration 8 cycles


Brady Settings Pacing Output 5.0 V @ 0.60 ms Mode SSI Sensitivity 1.0 mV Lower Rate Limit Lead Configuration 80 min⁻¹ Bipolar Maximum Tracking Rate -- min⁻¹ Pace Refractory 250 ms Sense

Daily Measurement

my measurement			
Date	Amplitude	Impedance	Threshold
	(mV)	(Ω)	(V)
09-JUL-2013	3.5	600	Off
08-JUL-2013	8.0	640	Off
07-JUL-2013	N.R.	610	Off
06-JUL-2013	N.R.	610	Off
05-JUL-2013	5.1	610	Off
04-JUL-2013	N.R.	620	Off
03-JUL-2013	N.R.	620	Off
27-JUN-2013	8.3	600	Off
20-JUN-2013	7.6	610	Off
13-JUN-2013	8.1	620	Off
06-JUN-2013	8.0	610	Off
30-MAY-2013	5.6	530	Off

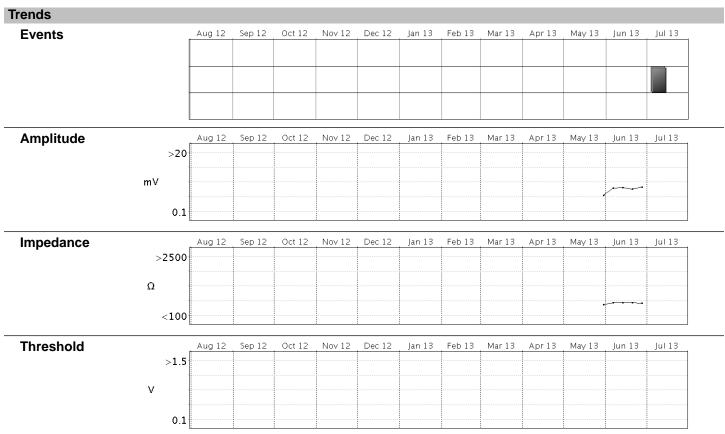

Tachy Counters		Since Last Reset		
Episode Counters				
Tachy Detections		5		
Brady Counters		Percentages	Since Last Reset	
Counters				
Paced		86 %	3.1M	
Sensed		14 %	498.5K	
Intrinsic Promotion				
Rate Hysteresis				
Searches			0	
Successful Searches			0	
Counters				
Interval Variation				
Variation 0 <= 10%			3.6M	
Variation 10 <= 20%			17.4K	
Variation 20 <= 30%			3.9K	
Variation > 30%			7.2K	
Histograms	Since Last Reset			

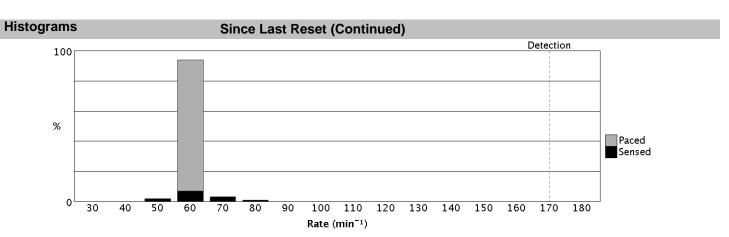

All Events Since Last Reset

09-JUL-2013 12:21	Tach at 190 min ⁻¹
09-JUL-2013 12:14	Tach at 280 min ⁻¹
09-JUL-2013 12:09	Tach at 240 min ⁻¹
09-JUL-2013 12:09	Tach at 260 min ⁻¹
09-JUL-2013 12:07	Tach at 220 min ⁻¹

2892 Software Version: 1.00 S508 Firmware Version: 0a © 2010

Boston Scientific Corporation or its affiliates. All rights reserved.

Page 3 of 3

Clinician Signature: