

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Overview		(1st interrog.)	
Mode		DDDR-ADIR	
Basic rate/Night rate [bpm]		60/50	
Sensor/Rate fading [bpm]		120/OFF	
Upper rate response [bpm]		130/2:1	
Mode switching [bpm]		160/DDIR	
Vp suppression		ON	
Dynamic AV delay [ms]		200/140	
MRI program		OFF	
Home Monitoring		OFF	
		A	V
Pulse amplitude [V]	⚡	3.6	⚡ 1.5
Pulse width [ms]		0.4	0.4
Capture control		ON	ON
Sensitivity [mV]		AUTO	AUTO
Refractory period/Blanking		Standard	
Sensing polarity	⚡	BIPL	⚡ BIPL
Pacing polarity	⚡	BIPL	⚡ BIPL
Calculated ERI		8 Y. 11 Mo.	

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Bradycardia	(1st interrog.)
---------------------------------	-----------------

Mode	DDDR-ADIR
Basic rate/Night rate [bpm]	60/50
Basic rate [bpm]	60
Night rate [bpm]	50
Night begins	22:00
Night ends	06:00
Hysteresis [bpm]	-5
Repetitive/Scan cycles	OFF
Atrial overdrive	OFF
Magnet response	AUTO
Sensor/Rate fading [bpm]	120/OFF
Max. activity rate [bpm]	120
Sensor gain	⚡ AUTO
Sensor threshold	Medium
Rate fading	OFF
Rate increase [bpm/cycle]	2
Rate decrease [bpm/cycle]	0.5
Upper rate response [bpm]	130/2:1
Upper tracking rate [bpm]	130
2:1 rate [bpm]	113
Atrial upper rate [bpm]	240
Mode switching [bpm]	160/DDIR
Mode switching	ON
Intervention rate [bpm]	160
Switch to	DDIR
Onset criterion [out of 8]	5
Resolution criterion [out of 8]	5
Change of basic rate [bpm]	+10
Rate stabilization during mode switching	OFF
2:1 Lock-in protection	ON
Vp suppression	ON
Pacing suppression [consecutive Vs]	6
Pacing support [out of 8 cycles]	3
MRI program	OFF

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Bradycardia	(1st interrog.)
---------------------------------	-----------------

AV delay		Individual	
After		Sense	Pace
at 60 bpm [ms]		155	200
at 80 bpm [ms]		125	170
at 100 bpm [ms]		115	160
at 120 bpm [ms]		105	150
at 140 bpm [ms]		95	140
Sense compensation [ms]			-45

AV hysteresis mode

AV hysteresis [ms]		
After		Pace
AV hysteresis at 60 bpm [ms]		450
AV hysteresis at 80 bpm [ms]		450
AV hysteresis at 100 bpm [ms]		450
AV hysteresis at 120 bpm [ms]		380
AV hysteresis at 140 bpm [ms]		300

AV repetitive/scan cycles

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Bradycardia		(1st interrog.)	
		A	V
Pulse amplitude [V]	⌀	3.6	⌀ 1.5
Pulse width [ms]		0.4	0.4
Capture control		ON	ON
Min. amplitude [V]		1.0	
Threshold test start [V]		3.0	3.0
Safety margin [V]		1.0	0.5
Search type	Time of day	Time of day	
Interval [h]			
Time of day		00:30	00:30
Status		OK	OK
Note A	Unstable rates present.		
Note V			
Last measured threshold [V]		2.6	1.0
Date		28/01/2021	05/02/2021
Time		00:35	10:01
Sensitivity [mV]		AUTO	AUTO
Sensing polarity	⌀	BIPL	⌀ BIPL
Pacing polarity	⌀	BIPL	⌀ BIPL
Refractory period/Blanking		Standard	
Atrial refractory period [ms]		AUTO	
Auto PVARP		ON	
PVARP [ms]	⌀	375	
PVARP after PVC [ms]		400	
2:1 rate [bpm]		113	
Ventricular refractory period [ms]		250	
Far-field protection after Vs [ms]		100	
Far-field protection after Vp [ms]		150	
Ven. blanking after Ap [ms]		30	
PMT protection		ON	
VA criterion [ms]		350	
Calculated ERI		8 Y. 11 Mo.	

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Home Monitoring	(1st interrog.)
-------------------------------------	-----------------

General settings	
PID	98
Home Monitoring	OFF
Time of transmission [hh:mm]	
High atrial rate	
Ongoing atrial episode	
High ventricular rate	
Event-based IEGM	

Last message	
Message type	---
Message created on	

Periodic IEGM for HM follow-up					
Cycle duration [days]					
	1st date	2nd date	3rd date	4th date	5th date
Date					

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Diagnostics		(1st interrog.)
Recording triggers		
High atrial rate		AT
High ventricular rate		ON
Patient trigger		OFF
Pre-trigger recording [%]		75
IEGM signal		Filtered
Statistics		
Start resting period [hh:mm]		02:00
Resting period duration [h]		4.0
AV delay adj. sensing test		ON
Thoracic impedance		ON
High atrial rate (HAR)		
HAR limit [bpm]		200
High ventricular rate (HVR)		
HVR limit [bpm]		180
HVR counter [Events]		8
Enable lead check		ON

Patient:	Follow-up on:
Device: Etrinsa 8 DR-T	05/02/2021
S/N: 68834692 (PID: 98)	11:17

Parameters - Patient	(1st interrog.)
-----------------------------	-----------------

ID	
Last name	
First name	
Date of birth	XX/XX/XXXX
Gender	XXX
Date of implant	11/08/2017
Hospital, City	
Physician	
Phone	
Symptom	XXX
ECG indication	
Etiology	
LVEF [%]	
NYHA	XXX
Remark	

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Parameters - Patient	(1st interrog.)
-----------------------------	-----------------

Leads

Lead model	
Manufacturer	
Serial number	
Type	
Implantation	XX/XX/XXXX
Channels	

Lead model	
Manufacturer	
Serial number	
Type	
Implantation	XX/XX/XXXX
Channels	

Lead model	
Manufacturer	
Serial number	
Type	
Implantation	XX/XX/XXXX
Channels	

Lead model	
Manufacturer	
Serial number	
Type	
Implantation	XX/XX/XXXX
Channels	

Patient:

Follow-up on:

Device: Etrinsa 8 DR-T

05/02/2021



S/N: 68834692 (PID: 98)

11:17


Diagnostics - Timing

(1st interrog.)


Rate histograms [% time]

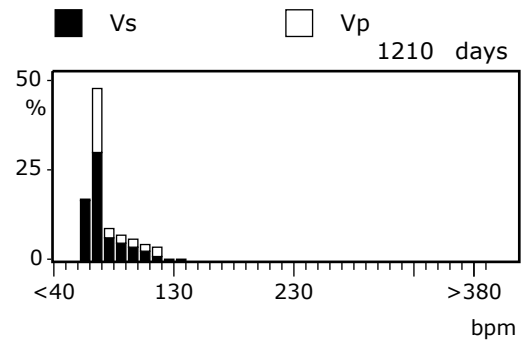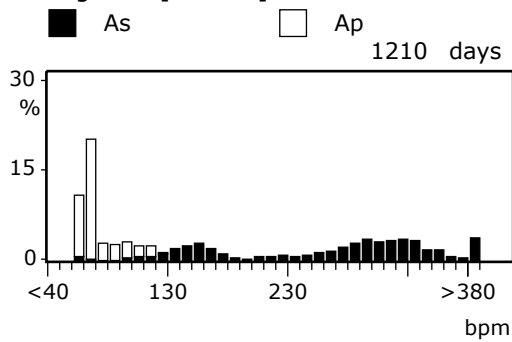


Ax-Vs histogram

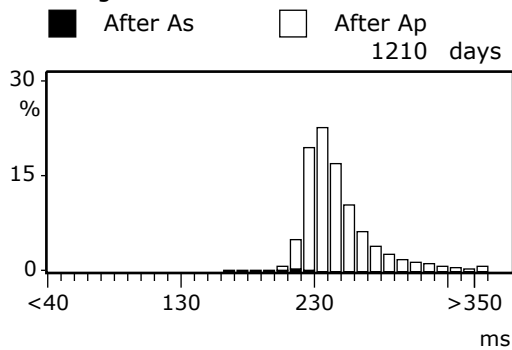

Event episodes and events

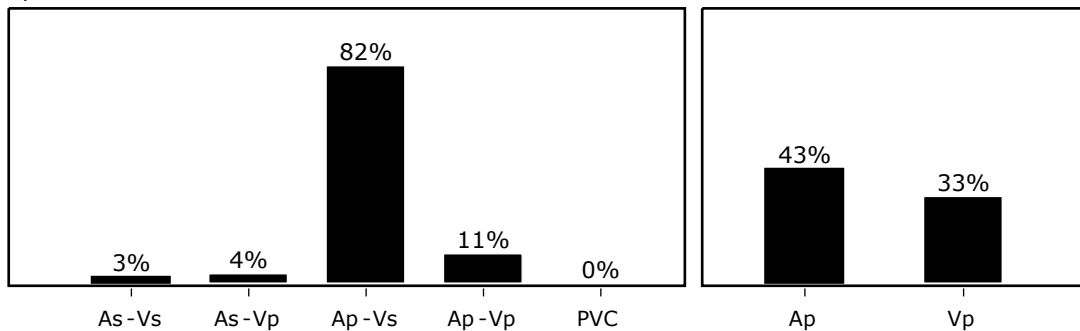

Pacing trends

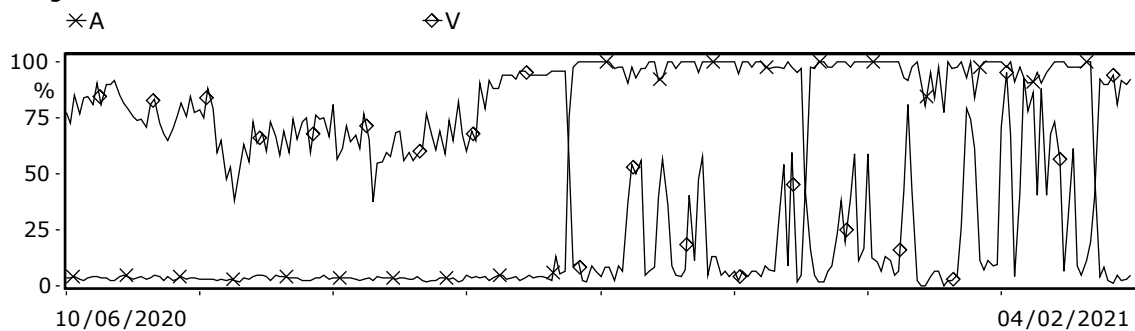

Start 13/10/2017

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	


Diagnostics - Arrhythmia A

(1st interrog.)


Atrial burden

Total number of episodes	12168
Atrial arrhythmia burden [%]	44
Longest atrial episode	0 days, 15:28


Time of occurrence

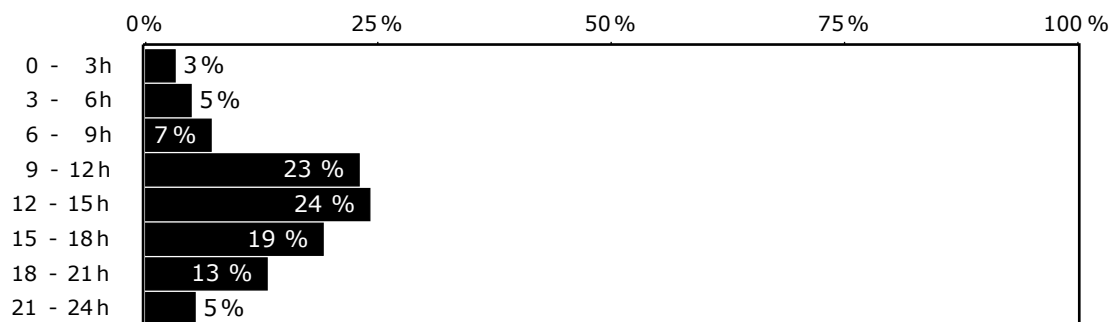

Tachy-episodes/24 h

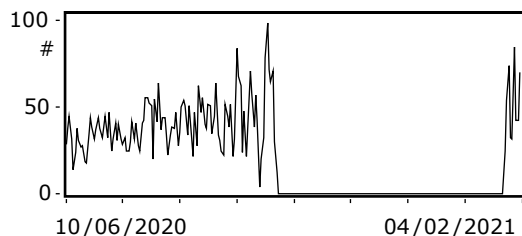

AT/AF burden

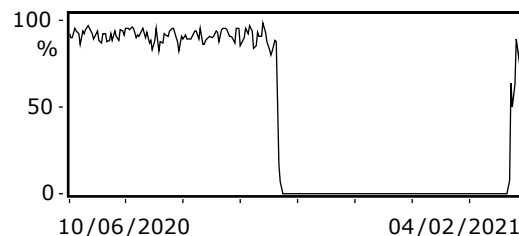

Mode switching

Mode switches [Count]	684063
Mode switching duration	287 days, 15:58

Start 13/10/2017

Patient:

Device: Etrinsa 8 DR-T

S/N: 68834692 (PID: 98)

Follow-up on:


05/02/2021

11:17


Diagnostics - HF monitor

(1st interrog.)


Mean rate

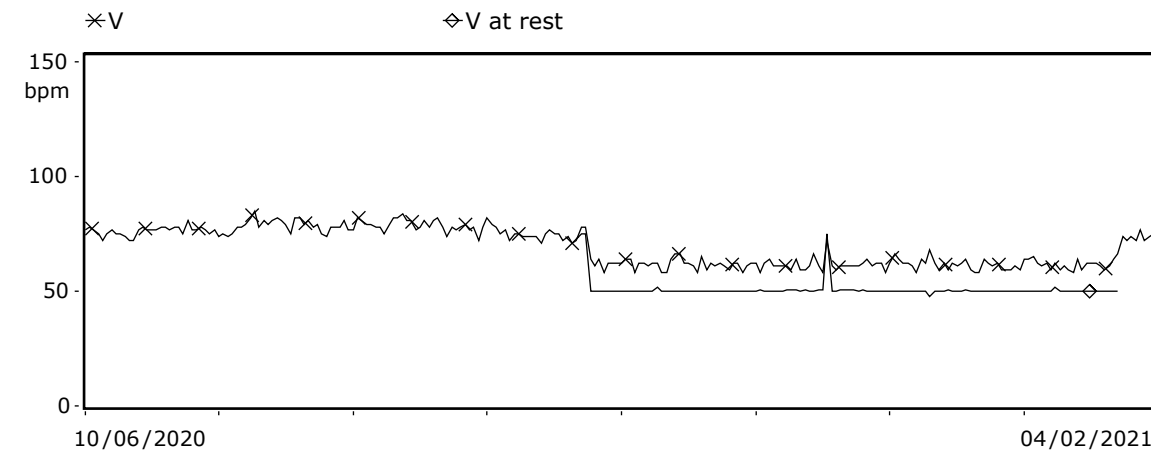

HR variability

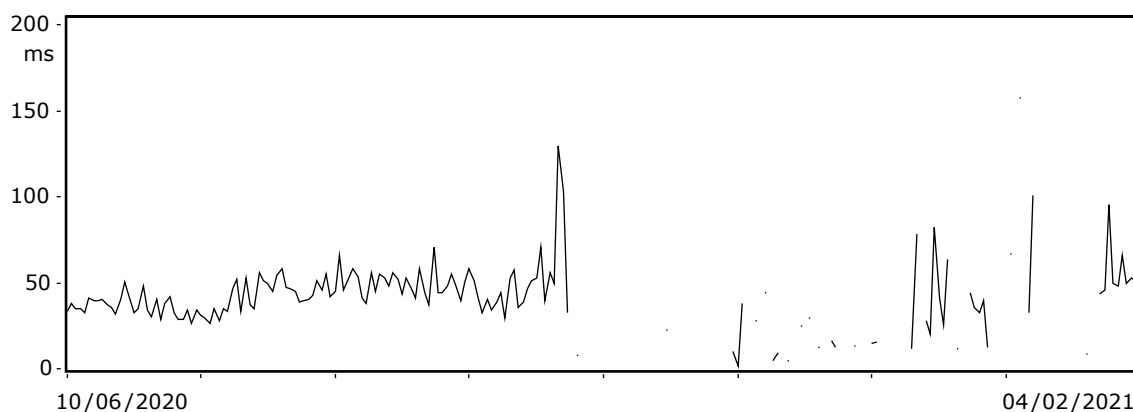

Patient activity

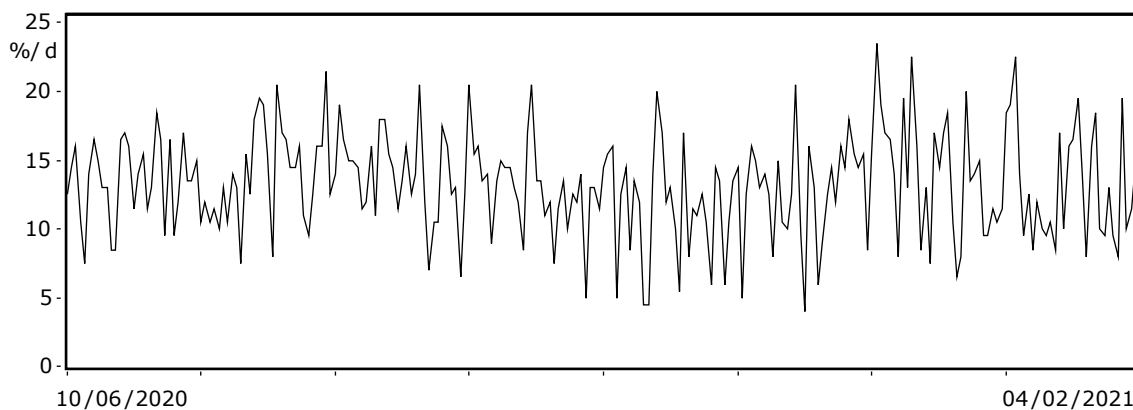


Start 13/10/2017

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Diagnostics - HF monitor	(1st interrog.)
---------------------------------	-----------------

Thoracic impedance

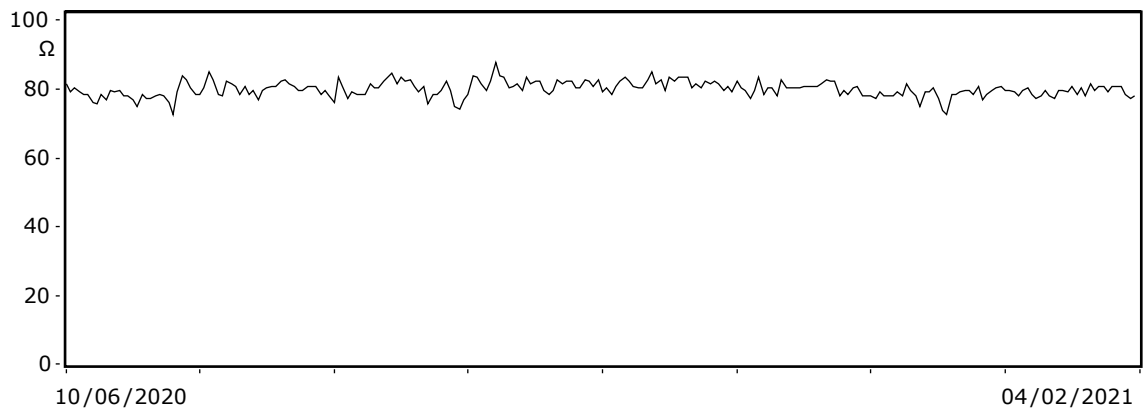

Start 13/10/2017

PSW 2201.A/1 / 5.0

12/26

.....
Patient:

.....
Follow-up on:

.....
Device: Etrinsa 8 DR-T

.....
05/02/2021


.....
S/N: 68834692 (PID: 98)

.....
11:17
.....


Diagnostics - 24 hours

(1st interrog.)

Rate

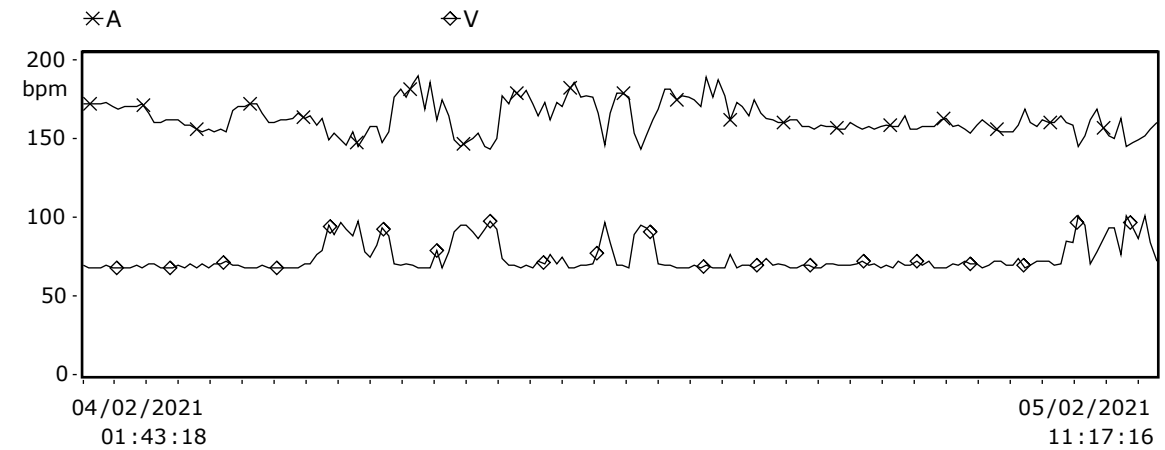

% paced

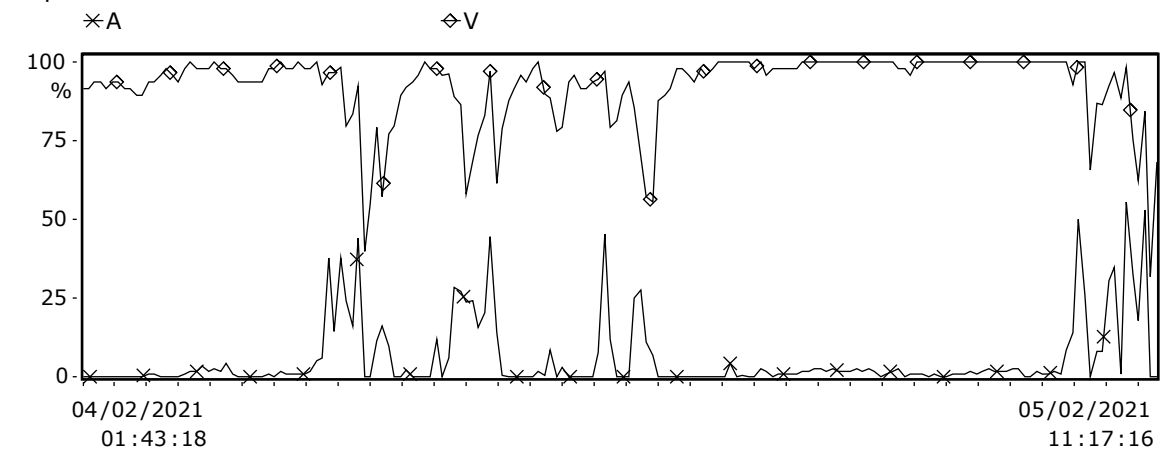

Start 13/10/2017

PSW 2201.A/1 / 5.0

13/26

Patient:

Follow-up on:

Device: Etrinsa 8 DR-T

05/02/2021



S/N: 68834692 (PID: 98)

11:17


Follow-up - Details

(1st interrog.)

Rate histograms [% time]

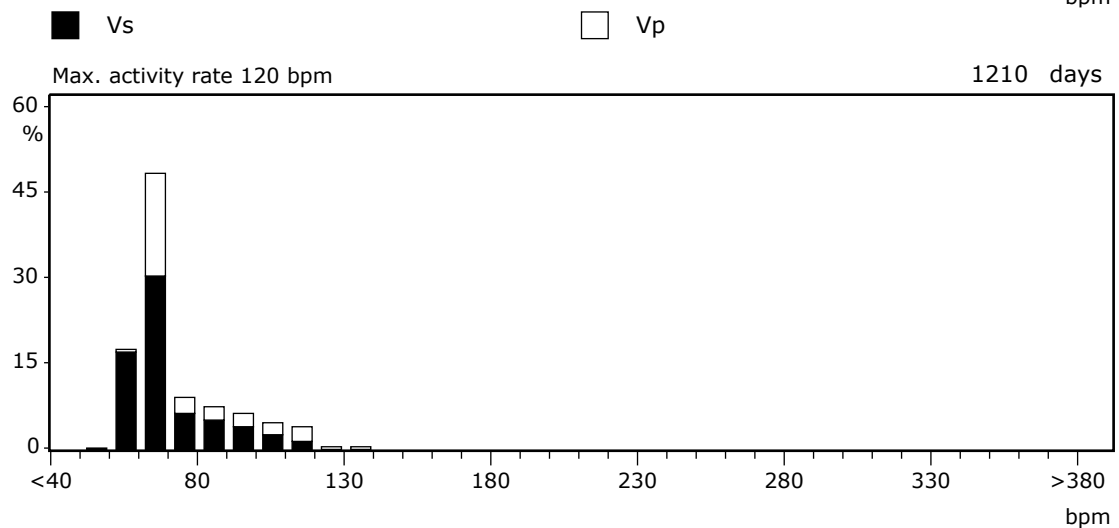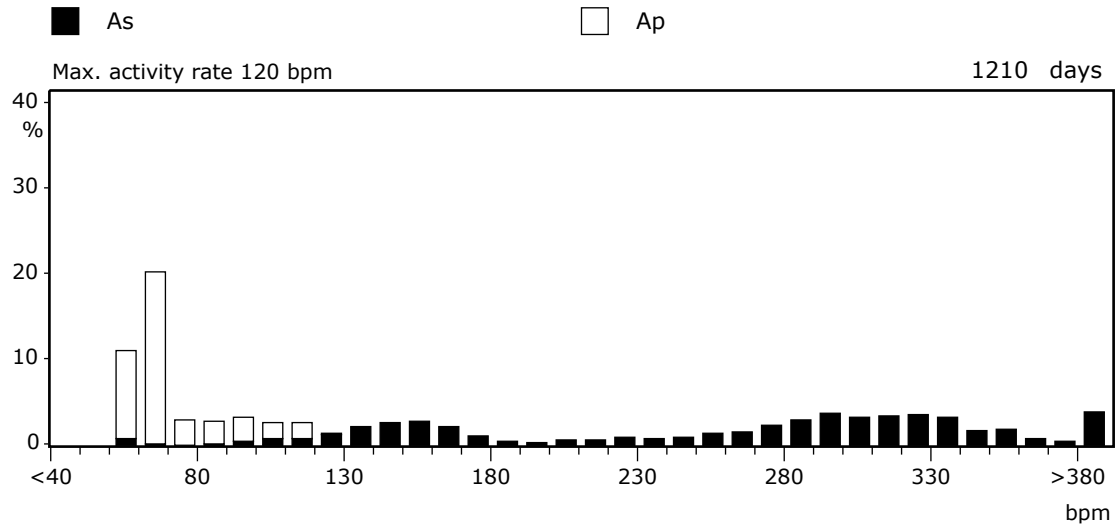

Event episodes

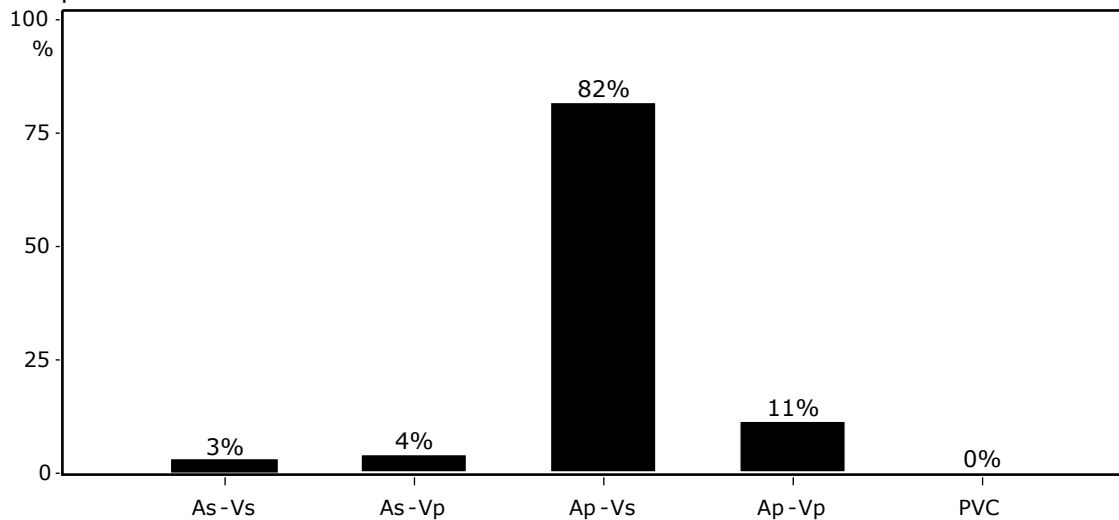


Start 13/10/2017

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	


Follow-up - Details	(1st interrog.)
----------------------------	-----------------

Long-term trends


Rate

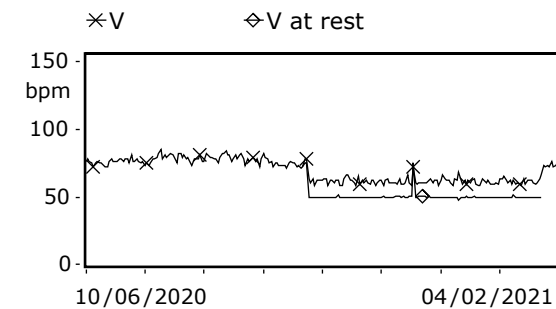


HR variability

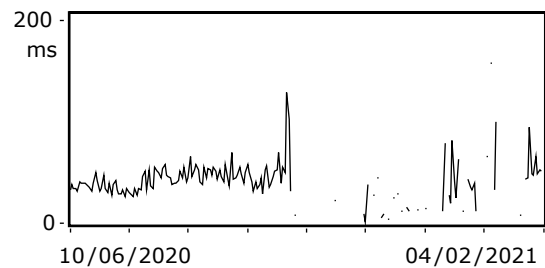

Patient activity

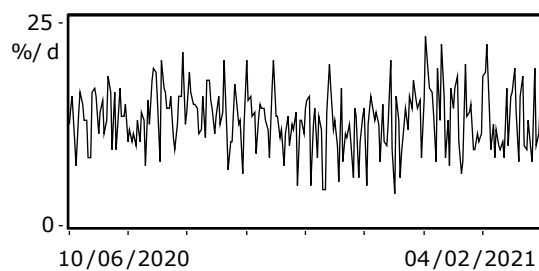

AT/AF burden

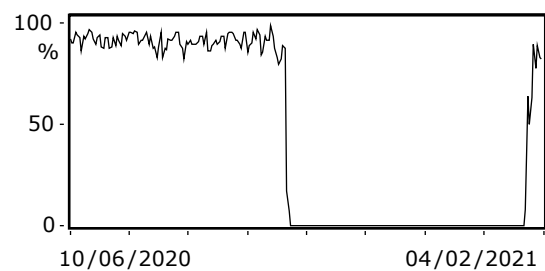

Patient:

Follow-up on:

Device: Etrinsa 8 DR-T

05/02/2021


S/N: 68834692 (PID: 98)

11:17


TrendView

(1st interrog.)


Mean rate

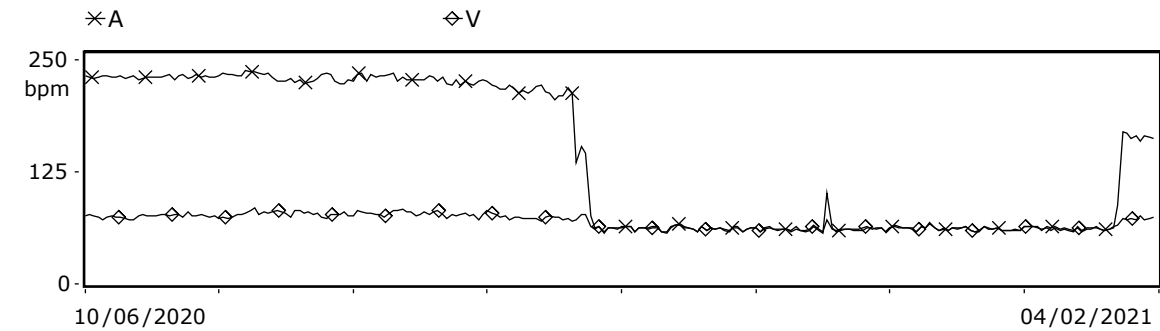


Pacing threshold

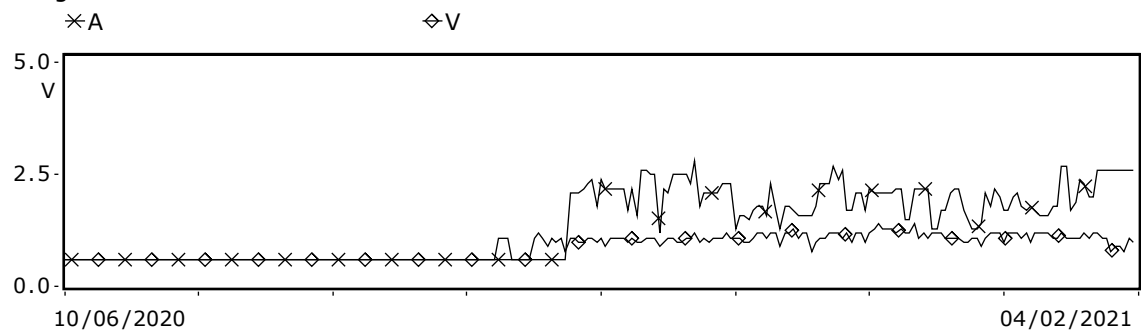

P/R wave amplitude

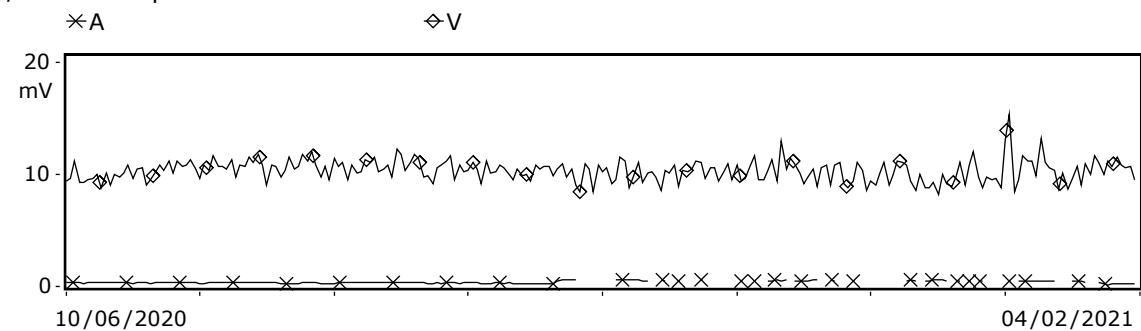

Lead impedance

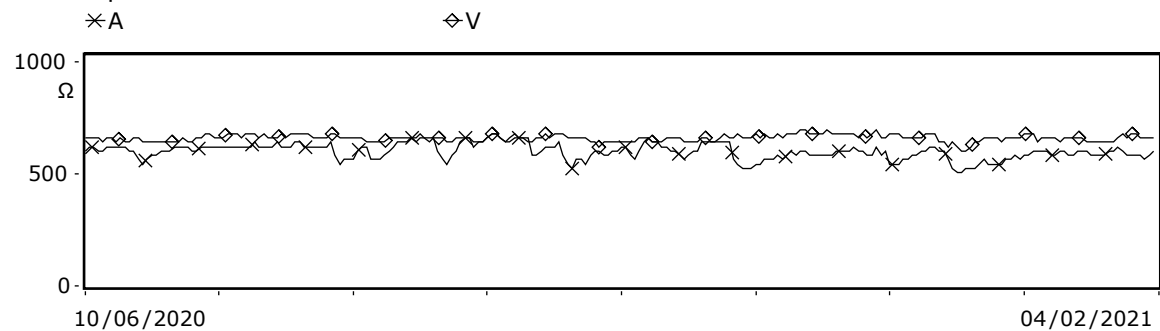




















Start 13/10/2017

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

TrendView		(1st interrog.)
Pacing in A [%]		43
Pacing in V [%]		33

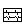Alerts and Episodes

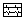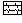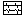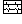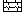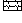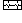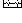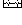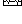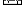Legend	
	Events
	IEGM recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording
	05/02/2021 AT recording

Patient:

Follow-up on:

Device: Etrinsa 8 DR-T

05/02/2021

S/N: 68834692 (PID: 98)

11:17

Recordings - Episodes

No.	Date	Duration	Trigger	MVR	PC
20	05/02/2021 12:28	00:31:12	AT	104	2
19	05/02/2021 12:10	00:05:34	AT	78	2
18	05/02/2021 12:06	00:00:22	AT	93	2
17	05/02/2021 11:55	00:01:34	AT	98	2
16	05/02/2021 11:32	00:13:04	AT	72	2
15	05/02/2021 11:22	00:02:26	AT	75	2
14	05/02/2021 11:17	00:03:46	AT	75	2
13	05/02/2021 11:07	00:02:40	AT	75	2
12	05/02/2021 11:03	00:02:24	AT	90	2
11	05/02/2021 11:00	00:01:24	AT	79	2
10	05/02/2021 10:36	00:20:58	AT	69	2
9	05/02/2021 10:32	00:00:38	AT	97	2
8	05/02/2021 10:29	00:00:52	AT	85	2
7	05/02/2021 10:04	00:08:22	AT	99	2
6	05/02/2021 09:37	00:20:30	AT	75	2
5	05/02/2021 09:09	00:27:38	AT	70	2
4	05/02/2021 08:28	00:32:50	AT	72	2
3	05/02/2021 08:16	00:11:10	AT	75	2
2	28/04/2020 17:43	* 15:28:46	AT	70	2
1	14/01/2019 18:40	01:01:48	AT	74	2

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:17	

Follow-up	(1st interrog.)
------------------	-----------------

Patient

Name	
Last follow-up	10/04/2019
Date of implant	11/08/2017

Device status

Program count	2
Home Monitoring	OFF

Mode	DDDR-ADIR			
Basic rate/UTR [bpm]	60/130			
		A		V
Pulse amplitude [V]	⌀	3.6	⌀	1.5
Pulse width [ms]		0.4		0.4
Expected ERI	8 Y. 11 Mo.			
Battery status	OK			
Remaining battery capacity [%]	75			

Test results

	A	V
Sensing amplitude [mV]	⌀ 0.4	⌀ 10.1
Sensing polarity	BIPL	BIPL
Threshold [V]		⌀ 1.0
Pacing polarity		BIPL
Pulse width [ms]		0.4
Lead impedance [Ω]	⌀ 624	⌀ 682
Pacing polarity	BIPL	BIPL

Diagnostics

Pacing in A/V [%]	43/33
Atrial arrhythmia burden [%]	44

.....
Patient:	Follow-up on:
.....
Device: Etrinsa 8 DR-T	05/02/2021
.....
S/N: 68834692 (PID: 98)	11:17
.....

Follow-up	(1st interrog.)
------------------	-----------------


Alerts and episodes

New episodes	20 recording(s)
--------------	-----------------

Patient: _____
Device: Etrinsa 8 DR-T
S/N: 68834692 (PID: 98)

Follow-up on: _____
05/02/2021
11:17

Recordings (05/02/2021 12:10 - AT (2))

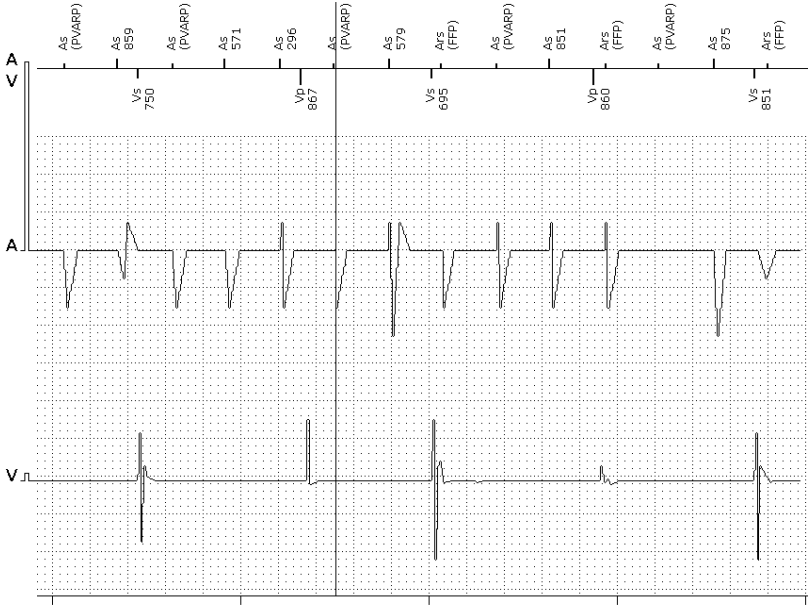


Follow-up on:

05/02/2021

11:17

Recordings

(05/02/2021 11:55 - AT (2))

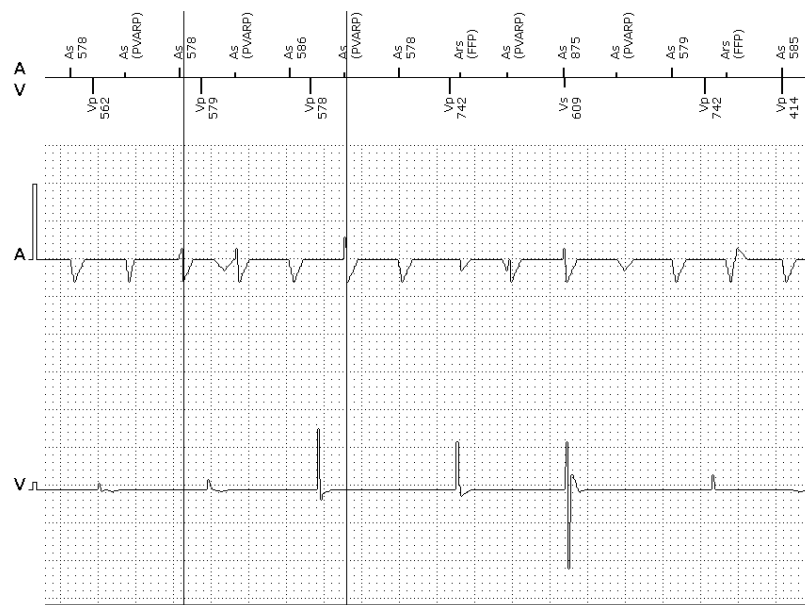

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:18	

Test - Sensing

Test program

Mode	P + R measurement DDI	
Basic rate [bpm]		30
AV delay [ms]		180
Upper tracking rate [bpm]		

	A	V
Sensitivity [mV]	AUTO	AUTO
Sensing polarity	BIPL	BIPL

Test settings

Print setup [mm/s]	OFF
--------------------	-----



Measured values

Mean rate [bpm]		199
	A	V
Min. amplitude [mV]	0.3	8.8
Mean amplitude [mV]	0.4	9.8
Max. amplitude [mV]	0.5	12.0

Follow-up on:

05/02/2021

11:18

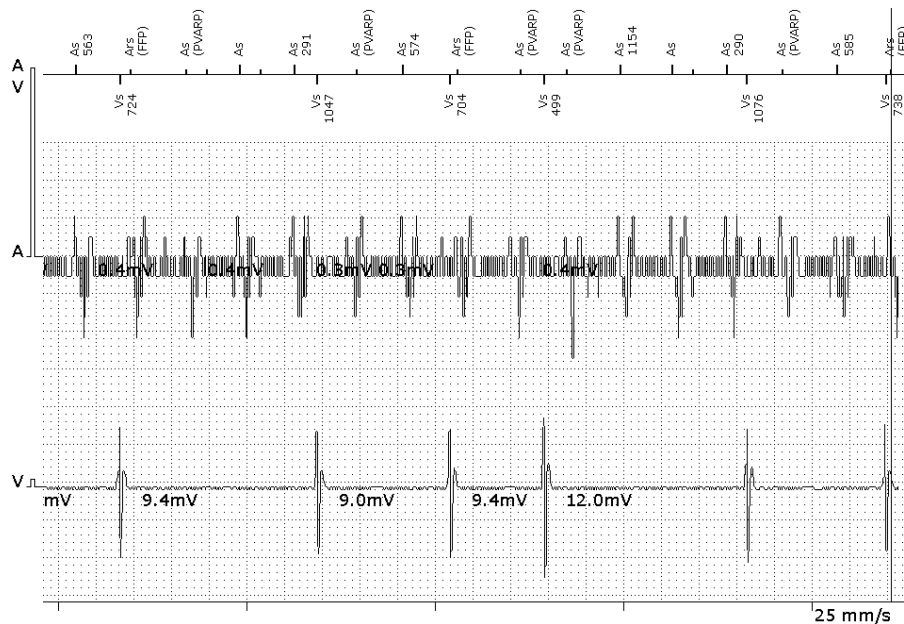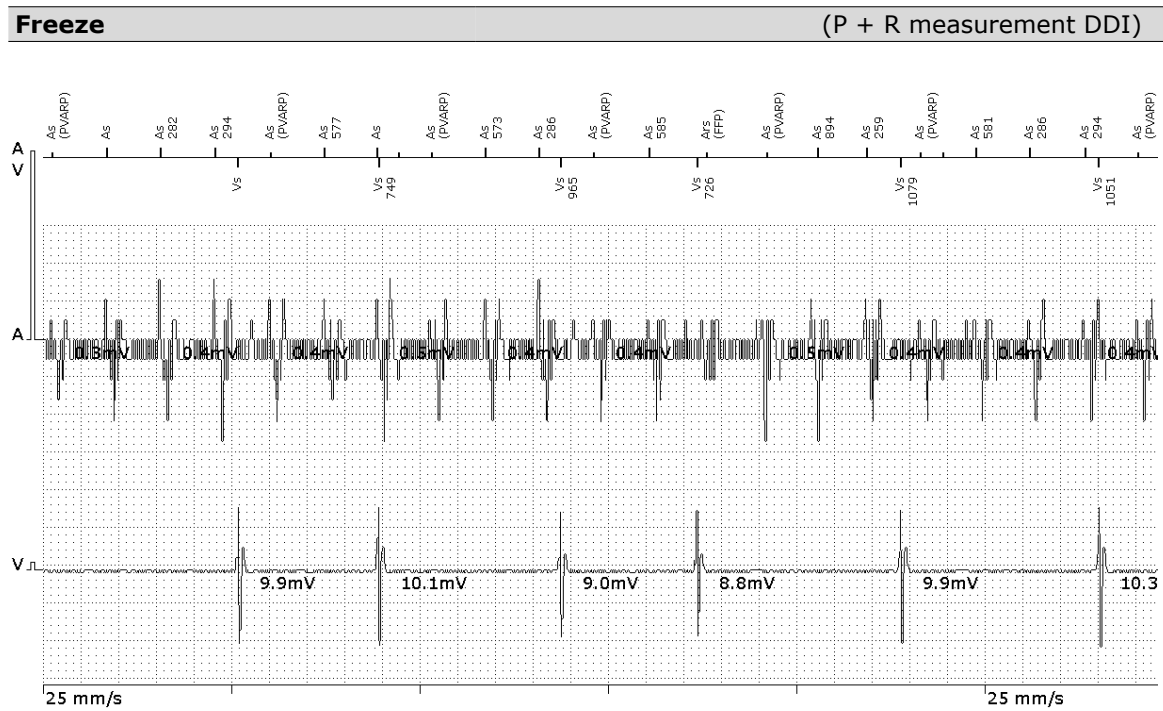

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:30	

Follow-up

Patient

Name	
Last follow-up	10/04/2019
Date of implant	11/08/2017

Device status

Program count	2
Home Monitoring	OFF

Mode	DDDR-ADIR		
Basic rate/UTR [bpm]	60/130		
	A	V	
Pulse amplitude [V]	⌀ 3.6	⌀ 1.5	
Pulse width [ms]	0.4	0.4	
Expected ERI	8 Y. 11 Mo.		
Battery status	OK		
Remaining battery capacity [%]	75		

Test results

	A	V	
Sensing amplitude [mV]	0.4	9.8	
Sensing polarity	BIPL	BIPL	
Threshold [V]		⌀ 1.0	
Pacing polarity		BIPL	
Pulse width [ms]		0.4	
Lead impedance [Ω]	⌀ 624	⌀ 682	
Pacing polarity	BIPL	BIPL	

Diagnostics

Pacing in A/V [%]	43/33
Atrial arrhythmia burden [%]	44

Patient:		Follow-up on:	
Device:	Etrinsa 8 DR-T	05/02/2021	
S/N:	68834692 (PID: 98)	11:30	

Follow-up

Alerts and episodes

New episodes	20 recording(s)
--------------	-----------------